Product Datasheet

SRC1 Antibody - BSA Free NB100-313

Unit Size: 100 ul

Store at 4C. Do not freeze.

www.novusbio.com



technical@novusbio.com

Publications: 1

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NB100-313

Updated 2/21/2025 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NB100-313



NB100-313


SRC1 Antibody - BSA Free

SIXCT Antibody - BOATTEE	
Product Information	
Unit Size	100 ul
Concentration	1.0 mg/ml
Storage	Store at 4C. Do not freeze.
Clonality	Polyclonal
Preservative	0.09% Sodium Azide
Isotype	IgG
Purity	Immunogen affinity purified
Buffer	Tris-Citrate/Phosphate (pH 7.0 - 8.0)
Product Description	
Host	Rabbit
Gene ID	8648
Gene Symbol	NCOA1
Species	Human
Immunogen	The immunogen this antibody was made to, maps to a region between residue 1400 and the C-terminus (residue 1441) of human Nuclear Receptor Coactivator 1 using the numbering given in entry NP_003734.3 (GeneID 8648).
Product Application Details	
Applications	Western Blot, Immunoprecipitation
Recommended Dilutions	Western Blot 1:1000-1:10000, Immunoprecipitation
Application Notes	NB 100-313 may be used for Western blot.




Images

Western Blot: SRC1 Antibody [NB100-313] - Detection of human and mouse SRC1 by western blot. Samples: Whole cell lysate from HEK293T and mouse NIH 3T3 cells prepared using NETN lysis buffer. Antibody: Affinity purified rabbit anti-SRC1 NB100-313 used for WB at 0.1 ug/ml. Detection: Chemiluminescence with an exposure time of 30 seconds.

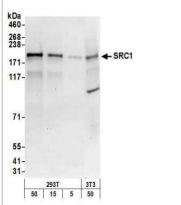

Immunoprecipitation: SRC1 Antibody [NB100-313] - Detection of human SRC1 by western blot of immunoprecipitates. Samples: Whole cell lysate (0.5 or 1.0 mg per IP reaction; 20% of IP loaded) from HEK293T cells prepared using NETN lysis buffer. Antibodies: Affinity purified rabbit anti-SRC1 antibodyNB100-313 (lot NB100-313-2) used for IP at 6 ug per reaction. SRC1 was also immunoprecipitated by a previous lot of this antibody (lot NB100-313-1). For blotting immunoprecipitated SRC1, NB100-313 was used at 0.4 ug/ml. Detection: Chemiluminescence with an exposure time of 30 seconds.

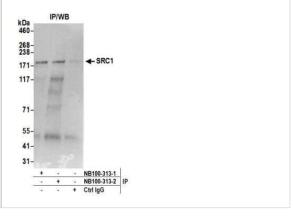

Publications

Heemers HV, Regan KM, Schmidt LJ, et al. Modulation of Coregulator Expression in Prostate Cancer Cells. Mol Endocrinol 2009-01-01 [PMID: 19164447]



Procedures

Serum protocol for SRC1 Antibody (NB100-313)

Nuclear Extract and Cytoplasmic Fraction Preparation protocol for SRC1 Antibody (NB100-313):

Nuclear Extract and Cytoplasmic Fraction Preparation

- 1. Nuclear extracts (NE) and cytoplasmic fractions (S100) were prepared by Dignam's method (Dignam, Lebovitz, and Roeder, Nucleic Acids Res. 11: 1475-1489. 1983).
- 2. 100 liters of HeLa cell culture were harvested and washed 3 times with cold PBS.
- 3. The packed-cell volume (PCV) was measured, and the cell pellet was gently resuspended with 5 PCVs of hypotonic buffer (10 mM HEPES-KOH [pH 8], 10 mM KCl, 1.5 mM MgCl2, 1 mM DTT, 0.2 mM PMSF).
- 4. Cells were incubated on ice for 10 minutes and then pelleted by centrifugation at 1,800xg for 10 minutes.
- 5. Hypotonic buffer was added to 2 PCVs, and cells were resuspended and then homogenized with 15 strokes using a pestle B in a Dounce glass homogenizer until the cells were more than 90% lysed, as determined by a light microscope.
- 6. The lysate was centrifuged at 20,000xg for 30 minutes at 4 degrees Celcius.
- 7. The supernatant was saved for S100 fraction, and the pellet was saved to measure the packed nuclear volume (PNV).
- 8. 0.4 ml of extraction buffer (20 mM HEPES-KOH [pH 8], 0.6 M KCl, 1.5 mM MgCl2, 0.2 mM EDTA, 25% [vol/vol] glycerol, 1 mM DTT, 0.2 mM PMSF) per ml of PNV was added.
- 9. Cell nuclei were homogenized with 10 strokes of pestle A in the homogenizer.
- 10. Suspension was stirred at 4 degrees Celcius for 30 minutes and centrifuged for 30 minutes at 20,000xg.
- 11. The supernatant (nuclear extract) was aliquotted for use.
- 12. The S100 fraction (resulting supernatant) was mixed with 0.11 volume of high-salt buffer (20 mM HEPES-KOH [pH 8], 1.2 M KCl, 1.5 mM MgCl2, 0.2 mM EDTA, 20% [vol/vol] glycerol, 1 mM DTT, 0.2 mM PMSF) and centrifuged at 100,000xg for 60 minutes at 4 degrees Celcius.
- 13. This supernatant was dialyzed for 2 hours at 4 degrees Celcius.
- 14. The sample was centrifuged for 30 minutes at 20,000xg and the supernatant (S100) was aliquotted for use.

Immunoprecipitation Antibody Characterization:

- 1. HeLa NE and S100 were diluted with 1 volume of RIPA buffer [150 mM NaCl, 1% NP-40, 0.5% DOC, 0.1% SDS, 50 mM Tris [pH 8]).
- 2. Cleared by spinning at 100,000 g for 20 minutes at 4 degrees Celcius.
- 3. 1 ml of supernatant (~10 mg total protein) was mixed with 20 ug of primary antibody (NB 100-313) and rotated overnight at 4 degrees Celcius.
- 4. Supernatant was mixed with 0.05 ml of protein A-sepharose beads (50% slurry) and rotated for 2 hours at 4 degrees Celcius.
- 5. Immunoprecipitates were washed 3 times with the 10% RIPA in PBS.
- 6. The washed beads were boiled with 0.04 ml of Laemmli buffer and subjected to SDS-PAGE (4-20% Tris-glycine gel).

Complex purification:

- 1. NE and S100 were cleared by spinning at 20,000 g for 30 minutes at 4 degrees Celcius.
- 2. 1.5 ml of supernatant (~15 mg total protein) was mixed with 20 ug of primary antibody (NB 100-313) and rotated for 4 hours at 4 degrees Celcius.
- 3. Sample and antibody mixture were centrifuged at 15,000 g for 20 minutes at 4 degrees Celcius.
- 4. Supernatant was mixed with 0.05 ml of protein A-sepharose beads (50% slurry) and rotated for 1 hour at 4 degrees Celcius.
- 5. Immunoprecipitates were washed 3 times with the NETN buffer (20 mM Tris-HCl [pH 8], 100 mM NaCl, 1 mM EDTA, 0.5% NP-40).
- 6. The washed beads were boiled with 0.04 ml of Laemmli buffer and subjected to SDS-PAGE (4-20% Tris-glycine gel).
- *If an insufficient amount of protein is purified for identification from 15 mg of extract, carry out the same procedure using 50-100 mg of extract to increase the amount of purified protein yield.





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112

USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com

Technical Support: nb-technical@bio-

techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Products Related to NB100-313

NB800-PC9 HeLa Nuclear Cell Lysate

HAF008 Goat anti-Rabbit IgG Secondary Antibody [HRP]

NB7160 Goat anti-Rabbit IgG (H+L) Secondary Antibody [HRP]

NBP2-24891 Rabbit IgG Isotype Control

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NB100-313

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

