Product Datasheet

GAPDH Antibody (6C5cc) - BSA Free NB600-502-0.2mg

Unit Size: 0.2 mg

Store at 4C. Do not freeze.

www.novusbio.com

technical@novusbio.com

Reviews: 3 Publications: 87

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NB600-502

Updated 2/21/2025 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NB600-502



NB600-502-0.2mg


GAPDH Antibody (6C5cc) - BSA Free

Product Information		
Unit Size	0.2 mg	
Concentration	Concentrations vary lot to lot. See vial label for concentration. If unlisted please contact technical services.	
Storage	Store at 4C. Do not freeze.	
Clonality	Monoclonal	
Clone	6C5cc	
Preservative	0.09% Sodium Azide	
Isotype	IgG1	
Purity	Protein A purified	
Buffer	PBS (pH 7.4)	
Target Molecular Weight	36 kDa	
Product Description		
Host	Mouse	
Gene ID	2597	
Gene Symbol	GAPDH	
Species	Human, Mouse, Rat, Porcine, Amphibian, Canine, Chinese Hamster, Feline, Fish, Rabbit, Bovine (Negative), Goat (Negative)	
Reactivity Notes	Please note that this antibody is reactive to Mouse and derived from the same host, Mouse. Additional Mouse on Mouse blocking steps may be required for IHC and ICC experiments. Please contact Technical Support for more information. Chinese Hamster reactivity reported in scientific literature (PMID: 26115091).	
Immunogen	Hybridoma clone has been derived from hybridization of Sp2/0 myeloma cells with spleen cells of Balb/c mice immunized with Rabbit GAPDH.	
Product Application Details		
Applications	Western Blot, Simple Western, ELISA, Immunoassay, Immunocytochemistry/ Immunofluorescence, Immunohistochemistry, Immunohistochemistry-Paraffin, Immunoprecipitation	
Recommended Dilutions	Western Blot 0.5-1 ug/ml, Simple Western 1:100, ELISA, Immunohistochemistry, Immunocytochemistry/ Immunofluorescence 1:10-1:2000, Immunoprecipitation, Immunohistochemistry-Paraffin, Immunoassay	
Application Notes	Use in IHC-P reported in scientific literature (PMID:34496231). Since GAPDH is expressed in all cells it is for example becoming the marker of choice for a loading control in Western Blotting. See <u>Simple Western Antibody Database</u> for Simple Western validation: tested in HeLa lysate (0.5 mg/ml); antibody dilution of 1:100; separated by size; detects a band at 42 kDa	



Images	
Western Blot: GAPDH Antibody (6C5) [NB600-502] - MUS81 promotes DNA synthesis during mitosis in BRCA2-deficient cells. Cell extracts were immunoblotted as indicated. SMC1 and GAPDH were used as loading controls. Image collected and cropped by CiteAb from the following publication (https://www.nature.com/doifinder/10.1038/ncomms15983), licensed under a CC-BY license.	e - DOX + DOX siPOLD3: + + + + + siMUS81: - + - + - + - + - + BRCA2 - 460 kDa POLD3 - 71 kDa SMC1 - 65 kDa MUS81 - 65 kDa GAPDH
Western Blot: GAPDH Antibody (6C5) [NB600-502] - Western blot analysis of GAPDH expression in 1) HeLa, 2) NTERA-2, 3) A-431, 4) HepG2, 5) MCF-7, 6) NIH 3T3, 7) PC-12 and 8) COS-7 whole cell lysates. Theoretical molecular weight: 36 kDa.	250> 150> 100> 75> 50> 37>
Western Blot: GAPDH Antibody (6C5) [NB600-502] - Western blot of FGF23 in the kidney. GAPDH antibody (6C5) served as an internal control. Each value represents the mean +/- SEM. Citation: Sugiura H, Matsushita A, Futaya M, Teraoka A, Akiyama K-i, Usui N, et al. (2018) Fibroblast growth factor 23 is upregulated in the kidney in a chronic kidney disease rat model. PLoS ONE 13(3): e0191706. https://doi.org/10.1371/journal.pone.0191706	(d) ad the week of
Simple Western: GAPDH Antibody (6C5) [NB600-502] - Simple Western lane view shows a specific band for GAPDH in 0.5 mg/ml of HeLa lysate. This experiment was performed under reducing conditions using the 12- 230 kDa separation system. Note: band observed higher than predicted 36 kDa molecular weight.	



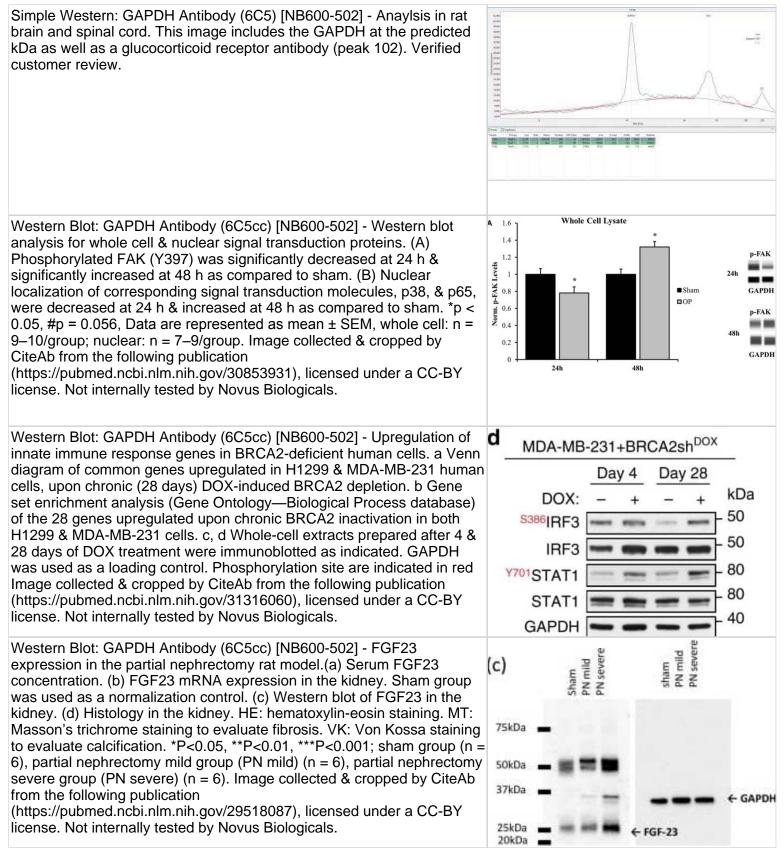



(d) Western Blot: GAPDH Antibody (6C5cc) [NB600-502] - FGF23 Wild type NP sham NP NX HP sham HP NX Wild type sham Nx sham Nx expression in hemi-nephrectomized rats fed a high-P diet.(a) Serum FGF23 concentration. (b) FGF23 mRNA expression in the bone. NP A A A H sham group was used as a normalization control. (c) FGF23 mRNA expression in the kidney. NP sham group was used as a normalization control. (d) Western blot of FGF23 in the kidney. GAPDH was used as an internal control. Each value shown represents the mean ± SEM; 50kDa *P<0.05, **P<0.01; NP sham group (n = 8), NP Nx group (n = 7), HP GAPDH 37kDa sham group (n = 7) & HP Nx group (n = 9). Image collected & cropped by CiteAb from the following publication FGF-23 25kDa (https://pubmed.ncbi.nlm.nih.gov/29518087), licensed under a CC-BY license. Not internally tested by Novus Biologicals. - DOX + DOX Western Blot: GAPDH Antibody (6C5cc) [NB600-502] - Time course of Day 18 18 innate immune response activation in BRCA2-deficient cells. a H1299 STING cells expressing a DOX-inducible BRCA2 shRNA were grown in the \$386IRF presence or absence of DOX for 28 days. Cells collected every 2 days IBES were subjected to quantitative RT-PCR analyses using primers specific STAT for indicated genes. mRNA levels were expressed relative to the gene STAT encoding β -actin & to day 2 (2- $\Delta\Delta$ CT). n = 3 independent experiments, GAPDH each performed in triplicates. Each dot represents a single replicate. b Whole-cell extracts prepared at indicated times after DOX addition were immunoblotted as shown. GAPDH was used as a loading control. Phosphorylation sites are indicated in red Image collected & cropped by CiteAb from the following publication (https://pubmed.ncbi.nlm.nih.gov/31316060), licensed under a CC-BY license. Not internally tested by Novus Biologicals. Western Blot: GAPDH Antibody (6C5cc) [NB600-502] - Decreased С proliferation rates in ALDH2 deficient human fibroblasts upon HR abrogationA-CHuman fibroblasts carrying either WT ALDH2 or the ALDH2 E487K homozygous ALDH2 mutation were transfected with control (SCR) or siRNAs against BRCA1 (A), BRCA2 (B), or RAD51 (C). Cells esiRAD5 were processed for proliferation assays 48 h after transfection. siRNAs were retransfected at an interval of 4 days. Graphs are representative of RAD5 three independent experiments, each performed with three technical replicas. Error bars represent SD of three technical replica values obtained from a single experiment. Cell extracts were prepared 48 h after GAPD transfection & immunoblotted as shown. SMC1 or GAPDH were used as loading controls. Image collected & cropped by CiteAb from the following publication (https://pubmed.ncbi.nlm.nih.gov/28729482), licensed under a CC-BY license. Not internally tested by Novus Biologicals. С Western Blot: GAPDH Antibody (6C5cc) [NB600-502] - Upregulation of H1299+BRCA2sh^{DOX} innate immune response genes in BRCA2-deficient human cells. a Venn diagram of common genes upregulated in H1299 & MDA-MB-231 human Day 4 Day 28 cells, upon chronic (28 days) DOX-induced BRCA2 depletion. b Gene DOX: kDa set enrichment analysis (Gene Ontology—Biological Process database) 50 S386IRF3 of the 28 genes upregulated upon chronic BRCA2 inactivation in both H1299 & MDA-MB-231 cells. c, d Whole-cell extracts prepared after 4 & 50 IRF3 28 days of DOX treatment were immunoblotted as indicated. GAPDH was used as a loading control. Phosphorylation site are indicated in red 80 Y701STAT1 Image collected & cropped by CiteAb from the following publication 80 (https://pubmed.ncbi.nlm.nih.gov/31316060), licensed under a CC-BY STAT1 license. Not internally tested by Novus Biologicals. 40 GAPDH




Page 4 of 8 v.20.1 Updated 2/21/2025

Western Blot: GAPDH Antibody (6C5cc) [NB600-502] - CtIP is not required for degradation of the single-stranded telomeric overhang generated by TRF2 inactivationImmortalized MEFs were infected with retroviruses expressing the indicated shRNAs, followed by selection with puromycin for 72 h. Cell extracts were prepared 48 h later & analysed by Western blotting as indicated. SMC1 & GAPDH were used as loading controls. *non-specific band.Mbol- & Alul-digested DNA from cells treated as in (A) was resolved by pulsed-field gel electrophoresis & probed with end-labelled (AACCCT)4 probe. Representative pulsed-field gel samples run under native & denatured conditions are shown.Quantification of the 3' overhang in cells treated as in (B). For each sample, the ss/total DNA ratios were expressed relative to the GFP shRNA-treated control. Error bars represent SD of two independent experiments. Image collected & cropped by CiteAb from the following publication (https://pubmed.ncbi.nlm.nih.gov/25582120), licensed under a CC-BY license. Not internally tested by Novus Biologicals.

Western Blot: GAPDH Antibody (6C5cc) [NB600-502] - ISG induction in HR-deficient human cells & tumors & olaparib impact on this process. a Upregulation of innate immune response genes in BRCA2-deleted ovarian tumors (n = 4) versus tumors with median BRCA2 mRNA expression (n = 145). b Upregulation of innate immune response genes in RAD51-deleted ovarian tumors (n = 11) versus tumors with median RAD51 mRNA expression (n = 145). Dots in graphs represent individual tumors. Middle line (white), median; box limits 25 & 75 percentiles; whiskers, minimum & maximum values. c H1299 cells expressing a DOX-inducible BRCA2 shRNA were grown in the presence or absence of DOX for 4 days. Then, olaparib (1 or 10 μ M) was added for 72 h, followed by quantitative RT-PCR analyses. Primers specific for the indicated genes were used. mRNA levels were expressed relative to the gene encoding β -actin & to untreated (-DOX) control cells (2- $\Delta\Delta$ CT). Error bars represent SD of n = 3 independent experiments. *, p < 0.05(unpaired two-tailed t test). d Whole-cell extracts from cells treated as in (c) were immunoblotted as shown. GAPDH & SMC1 were used as loading controls. Phosphorylation sites are indicated in red. e Cells treated as in (c) were pulse-labeled with EdU for 30 min. Frequency of cells in G1, S, & G2/M stages of the cell cycle were determined using FACS analyses of EdU incorporation & propidium iodide staining. Error bars represent SD of n = 3 independent experiments Image collected & cropped by CiteAb from the following publication (https://pubmed.ncbi.nlm.nih.gov/31316060), licensed under a CC-BY license. Not internally tested by Novus Biologicals.

Page 5 of 8 v.20.1 Updated 2/21/2025

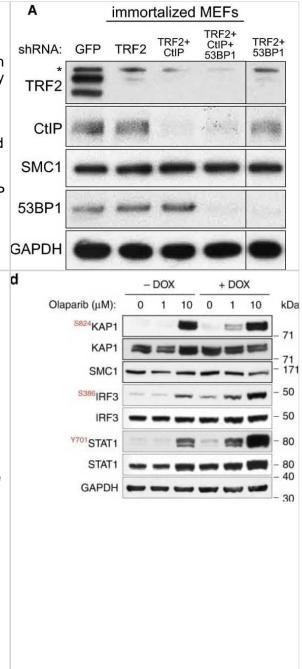


www.novusbio.com



Western Blot: GAPDH Antibody (6C5cc) [NB600-502] - BRCA1 & CtIP mediate Ku- & LIG4-independent telomere fusionsImmortalized Brca1F/- MEFs were infected with retroviruses expressing the indicated shRNAs and/or Cre recombinase, followed by selection with puromycin for 72 h. Mitotic chromosomes isolated 48 h later were fixed & stained with a Cy3-conjugated (CCCTAA)3-PNA probe. The frequency of end-toend chromosome-type fusions is represented as a percentage of fusions observed after TRF2 depletion. A minimum of 2,000 chromosomes were scored for each sample. Error bars represent SD of three independent experiments. P-values were calculated using an unpaired two-tailed ttest. NS, P > 0.05. Immortalized Brca1F/C61G MEFs were infected with retroviruses expressing the indicated shRNAs and/or Cre recombinase, followed by selection with puromycin for 72 h. The frequency of end-toend chromosome-type fusions was analysed as in (A).MEFs of the indicated genotypes were infected with retroviruses expressing TRF2 and/or CtIP shRNAs, followed by selection with puromycin for 72 h. Cell extracts were prepared 48 h later & analysed by Western blotting as indicated. GAPDH was used as a loading control. *non-specific band. Cells treated as in (C) & (E) were arrested in mitosis with colcemid, & mitotic chromosomes isolated 48 h later were fixed & stained with a Cy3conjugated (CCCTAA)3-PNA probe (D, F). The frequency of end-to-end chromosome-type fusions is represented as a percentage of fusions observed after TRF2 depletion. Error bars represent SD of two independent experiments. P-values were calculated using an unpaired two-tailed t-test. *P \leq 0.05. Image collected & cropped by CiteAb from the following publication (https://pubmed.ncbi.nlm.nih.gov/25582120), licensed under a CC-BY license. Not internally tested by Novus Biologicals.

Page 6 of 8 v.20.1 Updated 2/21/2025

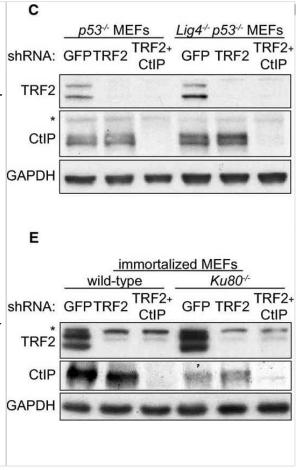



Publications

Park M, Baker W, Cambow D et al. Methamphetamine Enhances HIV-Induced Aberrant Proliferation of Neural Progenitor Cells via the FOXO3-Mediated Mechanism Molecular Neurobiology 2021-11-01 [PMID: 33983546] (Western Blot, Block/Neutralize)

Torices S, Teglas T, Naranjo O et al. Occludin Regulates HIV-1 Infection by Modulation of the Interferon Stimulated OAS Gene Family Molecular Neurobiology 2023-09-01 [PMID: 37209263] (Western Blot, Block/Neutralize)

Peyravian N, Sun E, Dikici E et al. Opioid Antagonist Nanodrugs Successfully Attenuate the Severity of Ischemic Stroke Molecular Pharmaceutics 2022-07-04 [PMID: 35506882]

Jones CE, Sharick JT, Colbert SE et al. Pten regulates collagen fibrillogenesis by fibroblasts through SPARC PloS one 2021-02-03 [PMID: 33534863]

Mendes VI, Bartholomeusz GA, Ayres M et al. Synthesis and cytotoxic activity of novel A-ring cleaved ursolic acid derivatives in human non-small cell lung cancer cells Eur J Med Chem 2016-07-22 [PMID: 27484517]

Dimitrios Kleidonas, Matthias Kirsch, Geoffroy Andrieux, Dietmar Pfeifer, Melanie Boerries, Andreas Vlachos Microglia modulate TNFα-mediated synaptic plasticity. Glia 2023-05-19 [PMID: 37208965]

Zakharova IO, Bayunova LV, Zorina II et al. Insulin and ?-Tocopherol Enhance the Protective Effect of Each Other on Brain Cortical Neurons under Oxidative Stress Conditions and in Rat Two-Vessel Forebrain Ischemia/Reperfusion Injury International Journal of Molecular Sciences 2021-10-29 [PMID: 34769198] (Western Blot)

van Senten JR, M
Iler TC, Moo EV et al. Use of CRISPR/Cas9-edited HEK293
cells reveals that both conventional and novel protein kinase C isozymes are involved in mGlu(5a) receptor internalization Journal of Biological Chemistry 2022-10-01 [PMID: 36087841] (Block/Neutralize)

Youn EK, Cho HM, Jung JK et al. Pathologic HDAC1/c-Myc signaling axis is responsible for angiotensinogen transcription and hypertension induced by high-fat diet Biomedicine & pharmacotherapy = Biomedecine & pharmacotherapie 2023-05-25 [PMID: 37244179] (WB, Mouse)

Torices S, Teglas T, Naranjo O et al. Occludin regulates HIV-1 infection by modulation of the interferon stimulated OAS gene family Research square 2023-01-30 [PMID: 36778388] (WB, Human)

Details:

Dilution used in WB 1:20,000

Zhao YM, Sun RS, Duan F et al. Intravitreal slow-release dexamethasone alleviates traumatic proliferative vitreoretinopathy by inhibiting persistent inflammation and M ller cell gliosis in rabbits International journal of ophthalmology 2023-01-18 [PMID: 36659954] (WB, Rabbit)

Details: Dilution used in WB 1:1000

Cao R, Chen P, Wang H et al. Intrafusal-fiber LRP4 for muscle spindle formation and maintenance in adult and aged animals Nature communications 2023-02-10 [PMID: 36765071] (WB, Mouse)

More publications at http://www.novusbio.com/NB600-502





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA Phone: 303.730.1950 Toll Free: 1.888.506.6887 Fax: 303.730.1966 nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402 canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449 Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com Orders: nb-customerservice@bio-techne.com General: novus@novusbio.com

Products Related to NB600-502-0.2mg

NBP1-42569	HepG2 Whole Cell Lysate
HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB720-B	Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin]
NBP1-97005-0.5mg	Mouse IgG1 Isotype Control (MG1)

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NB600-502

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

www.novusbio.com

