Product Datasheet

Profilin 1 Antibody - BSA Free NBP1-19344

Unit Size: 0.1 ml

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.



Publications: 1

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP1-19344

Updated 2/17/2025 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NBP1-19344




NBP1-19344

Profilin 1 Antibody - BSA Free

| Product Information | |
|-----------------------------|--|
| Unit Size | 0.1 ml |
| Concentration | 1.17 mg/ml |
| Storage | Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles. |
| Clonality | Polyclonal |
| Preservative | 0.1% Sodium Azide |
| Isotype | IgG |
| Purity | Immunogen affinity purified |
| Buffer | PBS and 30% Glycerol |
| Target Molecular Weight | 14 kDa |
| Product Description | |
| Host | Rabbit |
| Gene ID | 5216 |
| Gene Symbol | PFN1 |
| Species | Human, Mouse, Rat, Bovine |
| Immunogen | Synthetic peptide made to a C-terminal portion of human Profilin 1 (within residues 100-140). [Swiss-Prot# P07737] |
| Product Application Details | |
| Applications | Western Blot, Simple Western, Immunohistochemistry, Immunohistochemistry- Paraffin |
| Recommended Dilutions | Western Blot 0.5 ug/ml, Simple Western 1:200, Immunohistochemistry 1:100, Immunohistochemistry-Paraffin 1:100 |
| Application Notes | This Profilin 1 antibody is useful for Immunohistochemistry on formalin fixed paraffin embedded sections and Western blot, where a band is seen at approx. 14 kDa. |
| | In Simple Western only 10 - 15 uL of the recommended dilution is used per data point. See <u>Simple Western Antibody Database</u> for Simple Western validation: Tested in HeLa lysate 0.5 mg/mL, separated by Size, antibody dilution of 1:200, apparent MW was 18 kDa. Separated by Size-Wes, Sally Sue/Peggy Sue. |



Images

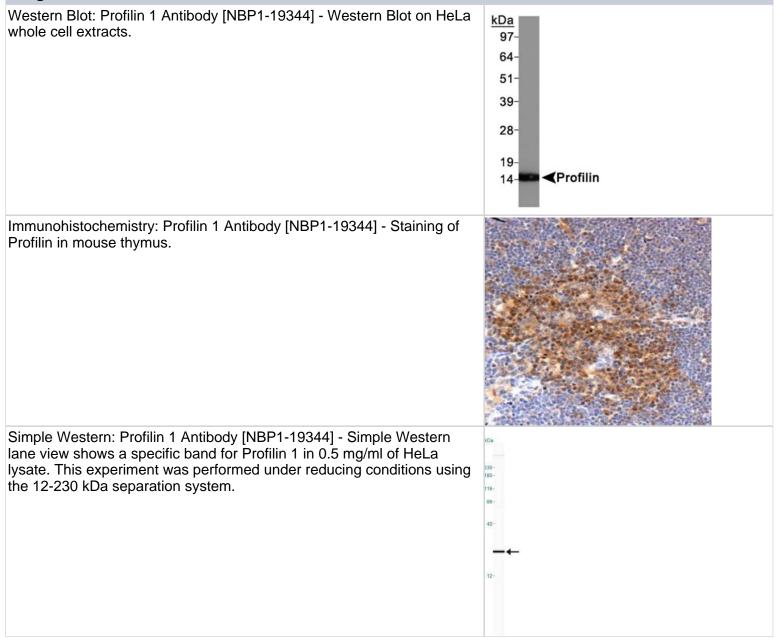

Publications

van Blitterswijk M, Baker MC, Bieniek KF et al. Profilin-1 mutations are rare in patients with amyotrophic lateral sclerosis and frontotemporal dementia. Amyotroph Lateral Scler Frontotemporal Degener 2013-05-02 [PMID: 23634771] (IHC-P, Human)



Procedures

Western Blot protocol for Profilin 1 Antibody (NBP1-19344) Profilin 1 Antibody:

Western Blot Protocol

1. Perform SDS-PAGE (4-12% MES) on samples to be analyzed, loading 30 ug of total protein per lane.

2. Transfer proteins to Nitrocellulose according to the instructions provided by the manufacturer of the transfer apparatus.

3. Rinse membrane with dH2O and then stain the blot using Ponceau S for 1-2 minutes to access the transfer of proteins onto the nitrocellulose membrane. Rinse the blot in water to remove excess stain and mark the lane locations and locations of molecular weight markers using a pencil.

4. Rinse the blot in TBS for approximately 5 minutes.

5. Block the membrane using 5% NFDM + 1% BSA in TBS + Tween, 1 hour at RT.

6. Rinse the membrane in dH2O and then wash the membrane in wash buffer [TBS + 0.1% Tween] 3 times for 10 minutes each.

7. Dilute the rabbit anti-Profilin primary antibody (NBP1-19344) in blocking buffer and incubate 1 hour at room temperature.

8. Rinse the membrane in dH2O and then wash the membrane in wash buffer [TBS + 0.1% Tween] 3 times for 10 minutes each.

9. Apply the diluted rabbit-IgG HRP-conjugated secondary antibody in blocking buffer (as per manufacturers instructions) and incubate 1 hour at room temperature.

10. Wash the blot in wash buffer [TBS + 0.1% Tween] 3 times for 10 minutes each (this step can be repeated as required to reduce background).

11. Apply the detection reagent of choice in accordance with the manufacturers instructions (Pierce ECL).

**Note: Tween-20 can be added to the blocking or antibody dilution buffer at a final concentration of 0.05-0.2%, provided it does not interfere with antibody-antigen binding.

Immunohistochemistry-Paraffin protocol for Profilin 1 Antibody (NBP1-19344)

Profilin 1 Antibody:

Immunohistochemistry-paraffin embedded sections

Antigen Unmasking

Bring slides to a boil in 10 mM sodium citrate buffer pH 6.0 then maintain at a sub-boiling temperature for 10 minutes. Cool slides on bench top for 30 minutes.

Staining

1. Wash sections in dH2O three times for 5 minutes each.

2. Wash section in wash buffer (1X PBS/0.1% Tween-20 (1X PBST)) for 5 minutes.

3. Block each section with 100-400 ul blocking solution (1X PBST, 5% goat serum) for 1 hour at room temperature.

4. Remove blocking solution and add 100-400 ul primary antibody diluted in 1X PBST, 5% goat serum to each section. Incubate overnight at 4C.

5. Remove antibody solution and wash sections in wash buffer three times for 5 minutes each.

6. Add 100-400 ul biotinylated secondary antibody, diluted in 1X PBST, 5% goat serum. Incubate 30 minutes at room temperature.

7. Remove secondary antibody solution and wash sections three times with wash buffer for 5 minutes each.

8. Add 100-400 ul Striptavidin-HRP reagent to each section and incubate for 30 minutes at room temperature.

9. Wash sections three times in wash buffer for 5 minutes each.

10. Add 100-400 ul DAB substrate to each section and monitor staining closely.

11. As soon as the sections develop, immerse slides in dH2O.

12. Counterstain sections in hematoxylin.

13. Wash sections in dH2O two times for 5 minutes each.

14. Dehydrate sections.

15. Mount coverslips.

www.novusbio.com





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA Phone: 303.730.1950 Toll Free: 1.888.506.6887 Fax: 303.730.1966 nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402 canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449 Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com Orders: nb-customerservice@bio-techne.com General: novus@novusbio.com

Products Related to NBP1-19344

| NB820-59177 | Human Brain Whole Tissue Lysate (Adult Whole Normal) |
|---------------|--|
| NBP1-19344PEP | Profilin 1 Antibody Blocking Peptide |
| HAF008 | Goat anti-Rabbit IgG Secondary Antibody [HRP] |
| NB7160 | Goat anti-Rabbit IgG (H+L) Secondary Antibody [HRP] |
| NBP2-24891 | Rabbit IgG Isotype Control |

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP1-19344

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

www.novusbio.com

