

# Product Datasheet

## GL7 (T and B Cell Activation Marker) Antibody (GL-7) [Biotin] NBP2-00361-50ug

Unit Size: 50 ug

Store at 4C. Do not freeze.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/NBP2-00361](http://www.novusbio.com/NBP2-00361)

Updated 5/15/2025 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/NBP2-00361](http://www.novusbio.com/reviews/destination/NBP2-00361)



**NBP2-00361-50ug**

GL7 (T and B Cell Activation Marker) Antibody (GL-7) [Biotin]


Product Information	
Unit Size	50 ug
Concentration	0.5 mg/ml
Storage	Store at 4C. Do not freeze.
Clonality	Monoclonal
Clone	GL-7
Preservative	0.09% Sodium Azide
Isotype	IgM
Conjugate	Biotin
Purity	Affinity purified
Buffer	PBS (pH 7.2)

Product Description	
Host	Rat
Species	Human, Mouse
Immunogen	The immunogen for this antibody is T- and B-cell activation antigen (GL7, Ly-77).

Product Application Details	
Applications	Flow Cytometry
Recommended Dilutions	Flow Cytometry 1:10-1:1000
Application Notes	The GL-7 (GL7) antibody has been tested by flow cytometric analysis of ConA-activated mouse splenocytes. This can be used at less than or equal to 0.5 ug per test. A test is defined as the amount (ug) of antibody that will stain a cell sample in a final volume of 100 uL. Cell number should be determined empirically but can range from $10^5$ to $10^8$ cells/test. The GL-7 antibody has also been reported for use in immunoprecipitation and immunohistochemical staining of frozen tissue sections. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

**Images**

Flow Cytometry: GL7 (T and B Cell Activation Marker) Antibody (GL-7) [Biotin] [NBP2-00361] - Staining of 3-day unstimulated (left) and 3-day ConA activated (right) BALB/c splenocytes with 0.5 ug of Rat IgM Isotype Control Biotin (open histogram) or 0.25ug of Anti-Human/Mouse GL7 (T and B Cell Activation Marker) Biotin (filled histogram) followed by Streptavidin PE. Total viable cells were used for analysis.

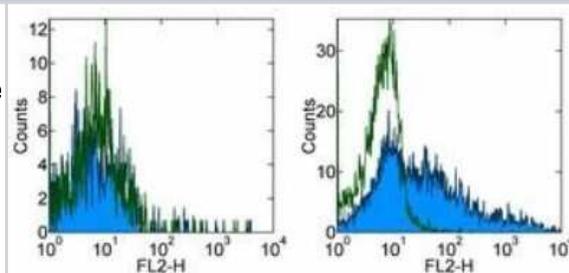



### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

www.novusbio.com  
Technical Support: nb-technical@bio-techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/NBP2-00361](http://www.novusbio.com/reviews/submit/NBP2-00361)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

