

Product Datasheet

NUT Antibody (2B6) - BSA Free NBP2-31147

Unit Size: 0.1 mg

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP2-31147

Updated 9/9/2025 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP2-31147



NBP2-31147


NUT Antibody (2B6) - BSA Free

Product Information	
Unit Size	0.1 mg
Concentration	1.0 mg/ml
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	2B6
Preservative	0.05% Sodium Azide
Isotype	IgG Kappa
Purity	Protein A purified
Buffer	PBS
Product Description	
Description	Novus Biologicals Mouse NUT Antibody (2B6) - BSA Free (NBP2-31147) is a monoclonal antibody validated for use in IHC, WB and ICC/IF. All Novus Biologicals antibodies are covered by our 100% guarantee.
Host	Mouse
Gene ID	256646
Gene Symbol	NUTM1
Species	Human, Mouse
Immunogen	Partial recombinant protein made to an internal portion of human NUT (between amino acids 400-700) [UniProt Q86Y26].
Product Application Details	
Applications	Western Blot, Immunohistochemistry-Paraffin, Immunocytochemistry/Immunofluorescence, Immunohistochemistry
Recommended Dilutions	Western Blot 2 ug/ml, Immunohistochemistry 2-5 ug/ml, Immunocytochemistry/Immunofluorescence 1:50 - 1:100, Immunohistochemistry-Paraffin 2-5 ug/ml




Images


Western Blot: NUT Antibody (46B4G9.2B6) [NBP2-31147] - WB analysis of partial recombinant NUT protein (expected molecular weight ~17 kDa) with NUT antibody (clone 46B4G9.2B6) at a concentration of 2 ug/ml.

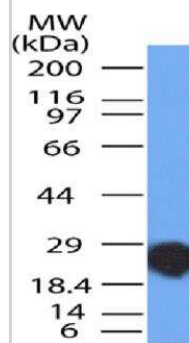

Immunocytochemistry/Immunofluorescence: NUT Antibody (46B4G9.2B6) [NBP2-31147] - NUT antibody was tested in A431 cells with Dylight 488 (green). Nuclei and alpha-tubulin were counterstained with DAPI (blue) and Dylight 550 (red). An antibody concentration of 0.01 ug/ml was used. Image objective 40x.

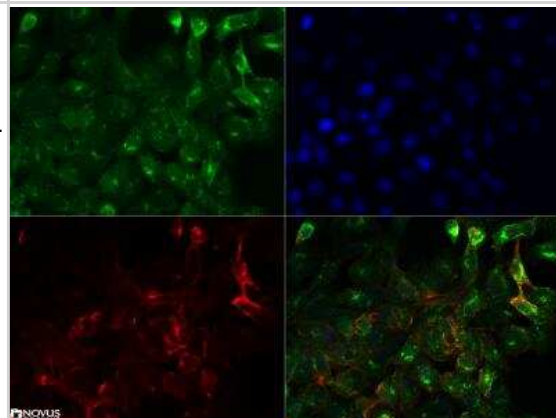

Immunohistochemistry-Paraffin: NUT Antibody (46B4G9.2B6) [NBP2-31147] - IHC-P analysis of NUT protein in normal human skin tissue section using NUT antibody (clone 46B4G9.2B6) at a concentration of 5 ug/ml. This antibody stained the nuclei and cytoplasm of different epidermal cell types such as granulosa cells, keratinocytes, Langerhans' cell and pigmented cells, and dermal blood vessels.

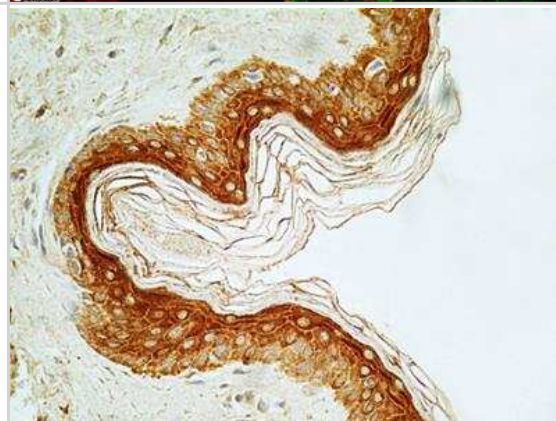

Procedures

Western Blot protocol for NUT Antibody (NBP2-31147)

NUT Antibody (2B6):

Western Blot Protocol

1. Perform SDS-PAGE on samples to be analyzed, loading 25 ug of total protein per lane.
 2. Transfer proteins to membrane according to the instructions provided by the manufacturer of the membrane and transfer apparatus.
 3. Stain according to standard Ponceau S procedure (or similar product) to assess transfer success, and mark molecular weight standards where appropriate.
 4. Rinse the blot.
 5. Block the membrane using standard blocking buffer for at least 1 hour.
 6. Wash the membrane in wash buffer three times for 10 minutes each.
 7. Dilute anti-NUT (2B6) primary antibody in blocking buffer and incubate 1 hour at room temperature.
 8. Wash the membrane in wash buffer three times for 10 minutes each.
 9. Apply the diluted HRP conjugated secondary antibody in blocking buffer (as per manufacturers instructions) and incubate 1 hour at room temperature.
 10. Wash the blot in wash buffer three times for 10 minutes each (this step can be repeated as required to reduce background).
 11. Apply the detection reagent of choice in accordance with the manufacturers instructions.
- Note: Tween-20 can be added to the blocking or antibody dilution buffer at a final concentration of 0.05-0.2%.

Immunocytochemistry/Immunofluorescence protocol for NUT Antibody (NBP2-31147)

NUT Antibody (2B6):

Immunocytochemistry Protocol

Culture cells to appropriate density in 35 mm culture dishes or 6-well plates.

1. Remove culture medium and add 10% formalin to the dish. Fix at room temperature for 30 minutes.
2. Remove the formalin and add ice cold methanol. Incubate for 5-10 minutes.
3. Remove methanol and add washing solution (i.e. PBS). Be sure to not let the specimen dry out. Wash three times for 10 minutes.
4. To block nonspecific antibody binding incubate in 10% normal goat serum from 1 hour to overnight at room temperature.
5. Add primary antibody at appropriate dilution and incubate at room temperature from 2 hours to overnight at room temperature.
6. Remove primary antibody and replace with washing solution. Wash three times for 10 minutes.
7. Add secondary antibody at appropriate dilution. Incubate for 1 hour at room temperature.
8. Remove antibody and replace with wash solution, then wash for 10 minutes. Add Hoechst 33258 to wash solution at 1:25,000 and incubate for 10 minutes. Wash a third time for 10 minutes.
9. Cells can be viewed directly after washing. The plates can also be stored in PBS containing Azide covered in Parafilm (TM). Cells can also be cover-slipped using Fluoromount, with appropriate sealing.

*The above information is only intended as a guide. The researcher should determine what protocol best meets their needs. Please follow safe laboratory procedures.

Immunohistochemistry-Paraffin protocol for NUT Antibody (NBP2-31147)**NUT Antibody (2B6):**

1. Deparaffinize the tissue sections by immersing the slides in Xylene with two changes for 10 min each. Sections should not get dried at any stage from this point.
2. Rehydrate the tissue sections by immersing the slides in decreasing grades of ethanol as follows:
 - a. Immerse in 100% ethanol with 2 changes for 5 minutes each
 - b. Immerse in 95% ethanol with 2 changes for 5 minutes each
 - c. Immerse in 90% ethanol for 5 minutes
 - d. Immerse in 70% ethanol for 5 minutes
 - e. Immerse in 50% ethanol for 5 minutes
 - f. Immerse in distilled water for 5 minutes
3. Antigen Retrieval (Microwave Method):
 - a. Immerse the slides in a microwave compatible tray containing 10 mM Sodium Citrate buffer (pH 6.0) with 0.05% Tween 20.
 - b. Boil the slides and maintain the sub-boiling temperature for 5 minutes in the microwave. Thereafter, take out the tray very carefully and cool it at room temperature (RT) for about 30 minutes.
 - c. Wash the slides 3 times, 3 minutes each by immersing them in TBST (Tris Buffered Saline having 0.05% Tween 20).
4. Quenching of Endogenous Peroxidase:
 - a. Incubate the slides in 3% hydrogen peroxide prepared in methanol for 15 minutes (at RT, in dark conditions).
 - b. Wash the slides in TBST 3 times, 3 minutes each.
5. Protein Blocking:
 - a. Incubate the sections with background sniper solution at RT for 15 minutes (Biocare Medicals, USA).
 - b. Wash the sections 3 times, 3 min each by immersing the slides in TBST.
6. Primary Antibody:
 - a. Dilute the primary antibody at 5ug/ml concentration using PBS as a diluent.
 - b. Incubate the sections with diluted primary antibody for 90 minutes at RT in a humidified chamber.
 - c. Thereafter, wash the slides 4 times, 5 minutes each with TBST.
7. Probe (Secondary Reagent):
 - a. Incubate with MACH 1 Mouse probe for 15 minutes at RT.
 - b. Incubate for 30 min at room temperature with HRP-Polymer (Biocare Medical, USA).
 - c. Wash the slides with TBST 4 times, 5 minutes each
8. Chromogen:
 - a. Mix 32ul of DAB Chromogen with 1 ml of DAB substrate buffer (Biocare Medical, USA).
 - b. Apply 200ul DAB mixture/section and incubate at RT in dark conditions (few seconds - 5 minutes).
 - c. As soon as an appropriate color develops, rinse the slides with deionized water (2-3 brief rinses).
9. Counter stain:
 - a. Counter stain with Hematoxylin for 30 seconds (Vector Labs, USA).
 - b. Wash in deionized water for 1-2 minutes to clear the extra stain.
 - c. Incubate the slides in bluing solution or Scott's water twice for 2 minutes each time.
10. Dehydrate the sections in increasing grades of alcohols:
 - a. 50% alcohol for 1 minute
 - b. 70% for 1 minute
 - c. 90% for 1 minute
 - d. 95% for 1 minute
 - e. 100% for 1 minute
 - f. Xylene with 2 changes for 2 minutes each
11. Mount with DPX mount and cover-slip glass (Fisher Scientific, USA), carefully not allowing any air bubbles to enter.

NOTE:- This protocol is provided as a reference tool only. Depending upon the type of tissues /tissue processing and reagents employed, the end user will need to optimize the final conditions for achieving an expected staining.



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Products Related to NBP2-31147

HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB720-B	Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin]
NBP1-84176PEP	NUT Recombinant Protein Antigen
NB100-56576	SARS Nucleocapsid Protein Antibody - BSA Free

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-31147

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

