## **Product Datasheet**

### Vimentin Antibody (VM452) - Azide and BSA Free NBP2-33060

Unit Size: 0.1 mg

Store at -20 to -80C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP2-33060

Updated 10/23/2024 v.20.1

# Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NBP2-33060



#### NBP2-33060


Vimentin Antibody (VM452) - Azide and BSA Free

Product Information	
Unit Size	0.1 mg
Concentration	1.0 mg/ml
Storage	Store at -20 to -80C. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	VM452
Preservative	No Preservative
Isotype	IgG1 Kappa
Purity	Protein A or G purified
Buffer	10 mM PBS
Target Molecular Weight	53.6 kDa
Product Description	
Description	<ul> <li>1.0 mg/ml of antibody purified from Bioreactor Concentrate by Protein A/G.</li> <li>Prepared in 10mM PBS WITHOUT BSA &amp; azide. Also available at 200 ug/ml WITH BSA &amp; azide (NBP2-29433).</li> <li>Antibody with azide - store at 2 to 8C. Antibody without azide - store at -20 to -80C.</li> </ul>
Host	Mouse
Gene ID	7431
Gene Symbol	VIM
Species	Human, Porcine, Bovine, Canine, Chicken, Feline, Goat, Mouse (Negative), Rat (Negative)
Reactivity Notes	Does not react with Mouse or Rat.
Marker	Mesenchymal Cell Marker
Specificity/Sensitivity	This monoclonal antibody reacts with a 58kDa protein identified as Vimentin. It shows no cross-reaction with other closely related intermediate filament proteins (IFP however, when used in panels with other antibodies, it is useful for the sub- classification of a given tumor. Expression of Vimentin, when used in conjunction with anti-keratin, is helpful when distinguishing melanomas from undifferentiated carcinomas and large cell lymphomas. All melanomas and Schwannomas react strongly with anti-Vimentin. It labels a variety of mesenchymal cells, including melanocytes, lymphocytes, endothelial cells, and fibroblasts. Non-reactivity of anti-Vimentin is often considered more useful than its positive reactivity, since there are a few tumors that do not contain Vimentin, e.g. hepatoma and seminoma. Anti-Vimentin is also useful as a tissue process control reagent.
Immunogen	Recombinant full-length human vimentin protein (Uniprot: P08670)
Product Application Details	
Applications	Western Blot, Flow Cytometry, Immunohistochemistry-Paraffin, Protein Array, Knockout Validated
Recommended Dilutions	Western Blot 0.5 - 2.0 ug/ml, Flow Cytometry, Immunohistochemistry-Paraffin, Protein Array, Knockout Validated
4	



#### Application Notes

Immunohistochemistry (Formalin-fixed): 1-2ug/ml for 30 minutes at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95C followed by cooling at RT for 20 minutes. Optimal dilution for a specific application should be determined.

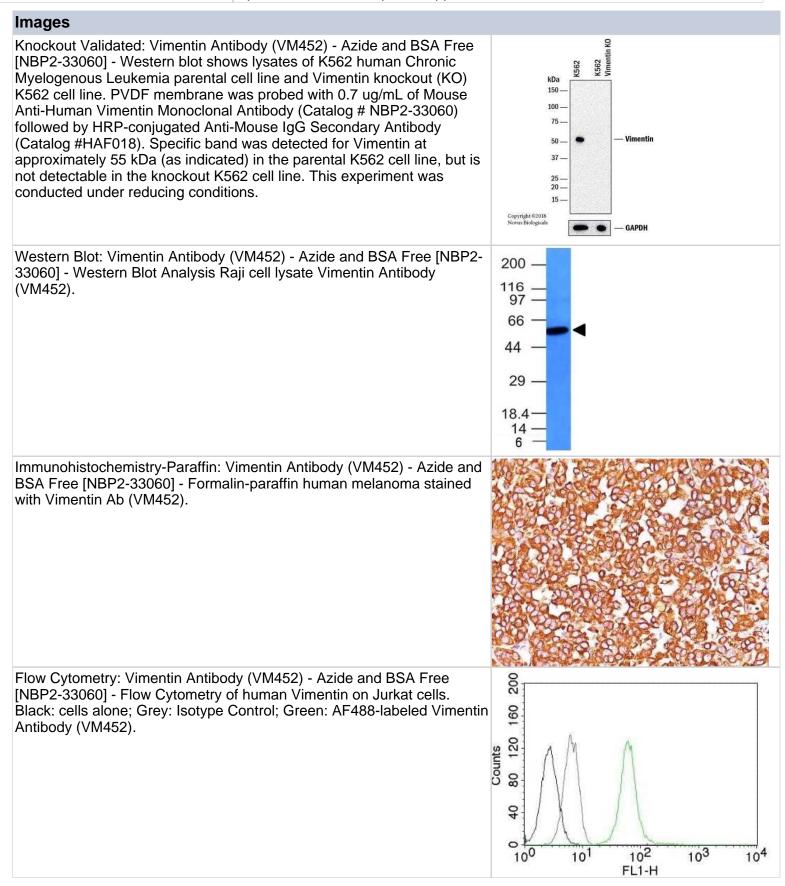

www.novusbio.com



Western Blot: Vimentin Antibody (VM452) - Azide and BSA Free [NBP2- 33060] - Western Blot Analysis of human A375 cell lysate using Vimentin Antibody (VM452).	
Protein Array: Vimentin Antibody (VM452) - Azide and BSA Free [NBP2- 33060] - Analysis of Protein Array containing more than 19,000 full- length human proteins using Vimentin Antibody (VM452). Z- and S- Score: The Z-score represents the strength of a signal that a monoclonal antibody (MAb) (in combination with a fluorescently-tagged anti-IgG secondary antibody) produces when binding to a particular protein on the HuProt(TM) array. Z-scores are described in units of standard deviations (SDs) above the mean value of all signals generated on that array. If targets on HuProt(TM) are arranged in descending order of the Z-score, the S-score is the difference (also in units of SDs) between the Z-score. S-score therefore represents the relative target specificity of a MAb to its intended target. A MAb is considered to specific to its intended target, if the MAb has an S-score of at least 2.5.	116         Postsin         Ruesto         8 score           87         2         CTES         1459         3 stat           90         3         MECh         163         4,09           90         1         10         20         30         40         50           90         1         10         20         30         40         50           Signal Rank (Top 50)         Signal Rank (Top 50)         50         50         50





#### Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA Phone: 303.730.1950 Toll Free: 1.888.506.6887 Fax: 303.730.1966 nb-customerservice@bio-techne.com

#### **Bio-Techne Canada**

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402 canada.inquires@bio-techne.com

#### **Bio-Techne Ltd**

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449 Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

#### **General Contact Information**

www.novusbio.com Technical Support: nb-technical@biotechne.com Orders: nb-customerservice@bio-techne.com General: novus@novusbio.com

#### Products Related to NBP2-33060

	Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin] Mouse IgG1 Kappa Isotype Control (P3.6.2.8.1)
NBP2-35139-100ug	Recombinant Human Vimentin Protein

#### Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-33060

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

www.novusbio.com

