## **Product Datasheet**

# HepG2 CoCl2 treated/Untreated Cell Lysate NBP2-36451

Unit Size: 2 Vials

Aliquot and store at -20C or -80C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP2-36451

Updated 10/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NBP2-36451



#### NBP2-36451

HepG2 CoCl2 treated/Untreated Cell Lysate


Product Information	
Unit Size	2 Vials
Concentration	1.0 mg/ml
Storage	Aliquot and store at -20C or -80C. Avoid freeze-thaw cycles.
Buffer	The protein lysate is prepared in 1x Laemmli sample buffer with BME. Boiling before loading is not necessary.
Decident December	

Product Description		
Species	Human	
Preparation Method	HepG2 cells are cultured using standard media condition until semi-confluent (70 -80%). The cells are then treated with or without 100 uM CoCl2 for 4 hours. The cells are then washed in PBS and lysed into 1X Laemmli sample buffer with betaME. The lysate is sonicated and tested by Western blot for reactivity to Hif-2a using a polyclonal anti-Hif-2a antibody (NB100-122). Tubulin reactivity of each lysate is shown as a loading control.	
Kit Components	0.1 ml Cobalt Chloride treated HepG2 whole cell lysate, 0.1 ml normoxic HepG2 whole cell lysate	
Lysate Type	Cell	
Lysate Tissue	Liver	
Lysate Subcellular Fraction	CoCl2 treated/Untreated	

-	
Product Application Details	
Applications	SDS-Page
Recommended Dilutions	SDS-Page
	HepG2 Whole cell lysate is provided as a Western blot positive control for Hif-2a research. The lysates are provided as a set of treated and untreated samples. Use 10 ul per lane for a standard mini-gel blot (approx. 1 mg/ml). Boiling before loading is not necessary.

## **Images**

Western Blot: HepG2 Hypoxic Simulated / Normoxic Tissue Lysate [NBP2-36451] - Western blot analysis of HIF-2 alpha antibody (NB100-122) in HepG2 cobalt chloride treated and normoxic lysates (NBP2-36451). Tubulin loading control is shown below.

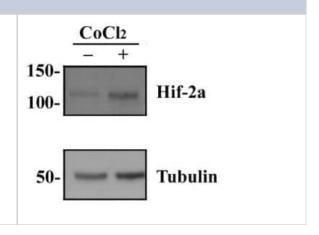





### **Novus Biologicals USA**

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

# Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

#### **Bio-Techne Canada**

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

#### **General Contact Information**

www.novusbio.com Technical Support: nb-technical@biotechne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

#### Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Lysates are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-36451

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

