Product Datasheet

TAG-72 Antibody (SPM536) - IHC-Prediluted NBP2-45176

Unit Size: 7 ml

Store at 4C.

www.novusbio.com

technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP2-45176

Updated 10/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NBP2-45176



NBP2-45176

TAG-72 Antibody (SPM536) - IHC-Prediluted


Product Information	
Unit Size	7 ml
Concentration	Please see the vial label for concentration. If unlisted please contact technical services.
Storage	Store at 4C.
Clonality	Monoclonal
Clone	SPM536
Preservative	0.05% Sodium Azide
Isotype	IgG1 Kappa
Purity	Protein A or G purified
Buffer	10 mM PBS with 0.05% BSA
Target Molecular Weight	220 kDa
Product Description	
Description	The prediluted antibody does not require any mixing, dilution, reconstitution, or titration; the antibody is ready-to-use and optimized for staining.
Host	Mouse
Species	Human, Rat, Bovine, Canine
Specificity/Sensitivity	Recognizes an oncofetal antigen of 220kDa, identified as a tumor-associated glycoprotein (TAG-72) with properties of a mucin. This monoclonal antibody defines the mucin-carried sialylated-Tn epitope. TAG-72 is usually expressed by adenocarcinomas, but is negative in mesotheliomas. Studies have reported that this antibody has 80% sensitivity and 93% specificity for pulmonary adenocarcinoma. Therefore, TAG-72 is a useful marker to distinguish between mesothelioma and adenocarcinoma. However, false positive reactions can occur so results must be interpreted with the utmost caution. This antibody may be useful in the differentiation of non-small cell carcinomas from small cell carcinomas of the lung. The combined use of anti-TAG-72 and anti-GCDFP-15 is valuable in the diagnosis of apocrine carcinoma.
Immunogen	Purified TAG-72 protein
Product Application Details	
Applications	Immunohistochemistry, Immunohistochemistry-Paraffin
Recommended Dilutions	Immunohistochemistry, Immunohistochemistry-Paraffin 0.5 - 1.0 ug/ml
Application Notes	Immunohistochemistry (Formalin-fixed): 1-2ug/ml for 30 minutes at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95C followed by cooling at RT for 20 minutes. Optimal dilution for a specific application should be determined.

www.novusbio.com



Images

Immunohistochemistry-Paraffin: TAG-72 Antibody (SPM536) - IHC-Prediluted [NBP2-45176] - Formalin-paraffin colon cancer stained with TAG-72 MAb (SPM536). Note cytoplasmic and cell surface staining.

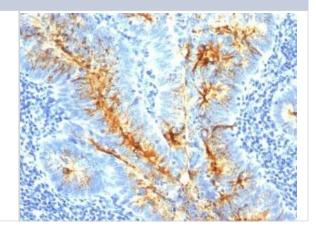





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA Phone: 303.730.1950 Toll Free: 1.888.506.6887 Fax: 303.730.1966 nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402 canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449 Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com Orders: nb-customerservice@bio-techne.com General: novus@novusbio.com

Products Related to NBP2-45176

HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB720-B	Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin]
NBP1-43319-0.5mg	Mouse IgG1 Kappa Isotype Control (P3.6.2.8.1)

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-45176

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

www.novusbio.com

