Product Datasheet

CD57 Antibody (HNK-1 + NK-1) - IHC-Prediluted NBP2-48350

Unit Size: 7 ml

Store at 4C.

www.novusbio.com

technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP2-48350

Updated 10/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NBP2-48350



NBP2-48350

CD57 Antibody (HNK-1 + NK-1) - IHC-Prediluted


7 ml		
Please see the vial label for concentration. If unlisted please contact technical services.		
Store at 4C.		
Monoclonal		
HNK-1 + NK-1		
0.05% Sodium Azide		
IgM Карра		
Protein L or PEG purified		
10 mM PBS with 0.05% BSA		
110 kDa		
Product Description		
The prediluted antibody does not require any mixing, dilution, reconstitution, or titration; the antibody is ready-to-use and optimized for staining.		
Mouse		
27087		
B3GAT1		
Human, Rat (Negative)		
Does not react with Rat.		
Natural Killer Cell Marker		
Anti-CD57 marks a subset of lymphocytes known as natural killer (NK) cells. Follicular center cell lymphomas often contain many NK cells within the neoplastic follicles. Anti-CD57 also stains neuroendocrine cells and their derived tumors, including carcinoid tumor and medulloblastoma. Anti-CD57 can also be useful in separating type B3 thymoma from thymic carcinoma when combined with a panel that includes antibodies against GLUT1, CD5, and CEA.		
Membrane antigen from HSB-2 cells		
Product Application Details		
Immunohistochemistry, Immunohistochemistry-Paraffin		
Immunohistochemistry, Immunohistochemistry-Paraffin		
Immunohistochemistry (Formalin-fixed): 2-4ug/ml for 30 minutes at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95C followed by cooling at RT for 20 minutes. Optimal dilution for a specific application should be determined.		

www.novusbio.com



Images

Immunohistochemistry-Paraffin: CD57 Antibody (HNK-1 + NK-1) - IHC-Prediluted [NBP2-48350] - Formalin-fixed, paraffin-embedded human tonsil stained with CD57 Monoclonal Antibody (HNK-1 + NK-1).

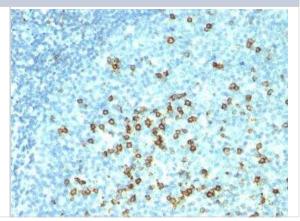

www.novusbio.com





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA Phone: 303.730.1950 Toll Free: 1.888.506.6887 Fax: 303.730.1966 nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402 canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449 Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com Orders: nb-customerservice@bio-techne.com General: novus@novusbio.com

Products Related to NBP2-48350

HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB720-B	Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin]
NBP1-96975-0.5mg	Mouse IgM Kappa Light Chain Isotype Control (MMK)
210-TA-005	TNF-alpha [Unconjugated]

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-48350

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

www.novusbio.com

