

Product Datasheet

Exosome Standards (COL01 cell line) NBP2-49845-200ug

Unit Size: 2 x 100ug Vials

Store at 4C. After reconstitution store at -70C.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP2-49845

Updated 10/23/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP2-49845



NBP2-49845-200ug

Exosome Standards (COLO1 cell line)

Product Information

Unit Size	2 x 100ug Vials
Concentration	Please see the protocols for proper use of this product. If no protocol is available, contact technical services for assistance.
Storage	Store at 4C. After reconstitution store at -70C.
Reconstitution Instructions	Add deionized water, 100 ul for Standard 100 ug and 30 ul for Standard 30 ug, to get a final concentration of 1 mg/mL. Resuspend exosomes pipetting the solution up and down 10-15 times, avoiding bubbles. Vortex the reconstituted standard for 60 seconds.
Buffer	Lyophilized from cell culture media

Product Description


Description	Highly pure, lyophilized exosome standards with superior stability, optimal for multiple applications including: Assay calibration, Spike-in control for exosome quantification, Protein marker analysis for different techniques such as Western Blot and Flow Cytometry, Extraction and analysis of exosomal RNA and DNA. Quantity per vial of 30 ug size (number of particles in 30 ug: $> 1 \times 10^8$). Quantity per vial of 100 ug size (number of particles in 100 ug: $> 1 \times 10^{10}$).
Preparation Method	Isolation involves Tangential flow filtration combined with Size Exclusion Chromatography. Exosomes (small EVs) are quantified and validated for protein content and particle number by Nanoparticle Tracking Analysis as well as for common tetraspanins marker validation. Lyophilization does not alter stability of exosome proteins and nucleic acids.

Product Application Details

Applications	Western Blot, ELISA, Electron Microscopy, Flow Cytometry, Nucleic Acid Extraction
Recommended Dilutions	Western Blot, Flow Cytometry, ELISA, Electron Microscopy, Nucleic Acid Extraction

Images

Western Blot: Exosome Standards (COLO1 cell line) [NBP2-49845] - Marker detection by Western Blotting. Reconstituted Exosomes can be directly lysed in Laemmli buffer, then loaded on the Electrophoresis gel. Recommended quantity: 10-20 ug per line.

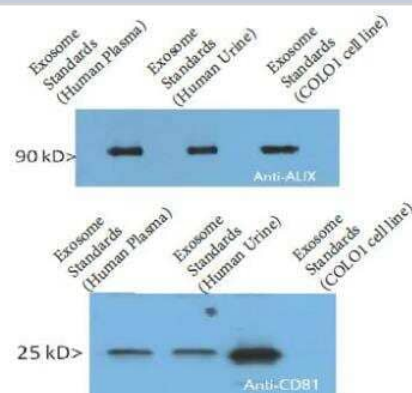



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-
techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Support products are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-49845

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

