

Product Datasheet

Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00) NBP2-67295

Unit Size: 100 ul

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Publications: 8

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP2-67295

Updated 10/23/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP2-67295



NBP2-67295

Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00)

Product Information

Unit Size	100 ul
Concentration	1 mg/ml
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	SA40-00
Preservative	0.05% Sodium Azide
Isotype	IgG
Purity	Protein A purified
Buffer	TBS (pH7.4), 0.05% BSA, 40% Glycerol

Product Description


Host	Rabbit
Gene ID	249
Gene Symbol	ALPL
Species	Human, Mouse, Rat
Immunogen	Synthetic peptide within human Alkaline Phosphatase, Tissue Non-Specific aa 18-50. (SwissProt: P05186 Human; SwissProt: P09242 Mouse; SwissProt: P08289 Rat)

Product Application Details


Applications	Western Blot, Flow Cytometry, Immunocytochemistry/ Immunofluorescence, Immunohistochemistry, Immunohistochemistry-Paraffin, Immunoprecipitation
Recommended Dilutions	Western Blot 1:5,000, Flow Cytometry 1:50, Immunohistochemistry, Immunocytochemistry/ Immunofluorescence 1:50 , Immunoprecipitation Use at an assay dependent concentration, Immunohistochemistry-Paraffin 1:50-1:200

Images


Immunocytochemistry/Immunofluorescence: Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00) [NBP2-67295] - Staining ALP in SW480 cells (green). Cells were fixed in paraformaldehyde, permeabilised with 0.25% Triton X100/PBS.

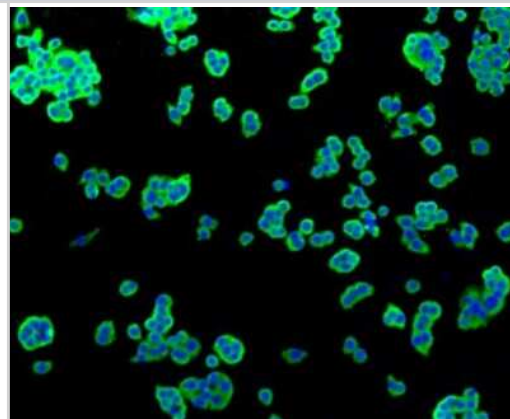


Flow Cytometry: Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00) [NBP2-67295] - Analysis of Hela cells with ALP antibody at 1/50 dilution (blue) compared with an unlabelled control (cells without incubation with primary antibody; red). Alexa Fluor 488-conjugated goat anti rabbit IgG was used as the secondary antibody

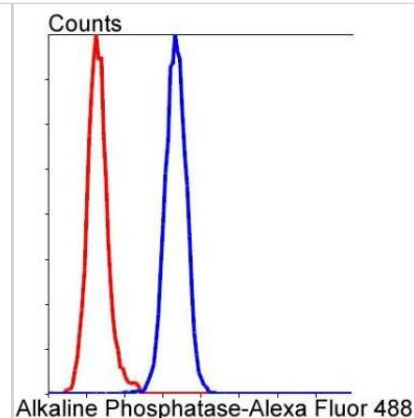


Immunocytochemistry/Immunofluorescence: Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00) [NBP2-67295] - Staining ALP in HepG2 cells (green). The nuclear counter stain is DAPI (blue). Cells were fixed in paraformaldehyde, permeabilised with 0.25% Triton X100/PBS.

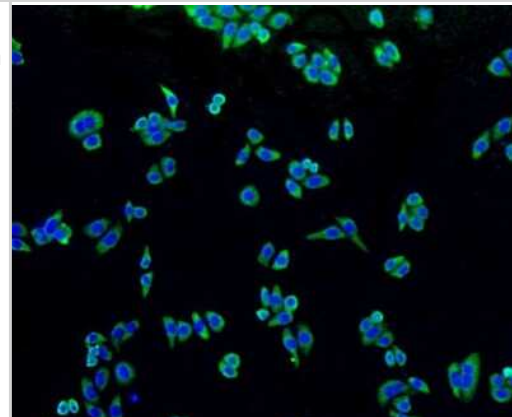


Immunohistochemistry-Paraffin: Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00) [NBP2-67295] - Immunohistochemical analysis of paraffin-embedded human kidney tissue using anti-Alkaline Phosphatase, Tissue Non-Specific antibody. The section was pre-treated using heat mediated antigen retrieval with Tris-EDTA buffer (pH 8.0-8.4) for 20 minutes. The t

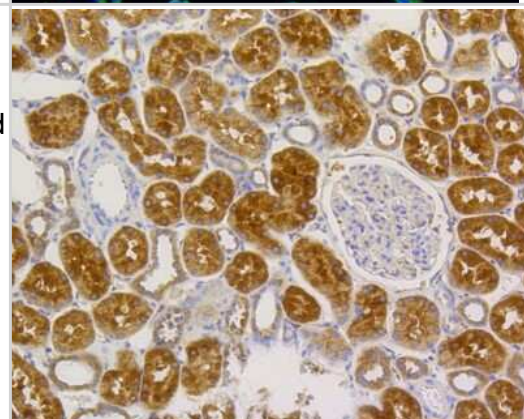


Immunohistochemistry-Paraffin: Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00) [NBP2-67295] - Immunohistochemical analysis of paraffin-embedded human liver carcinoma tissue using anti-Alkaline Phosphatase, Tissue Non-Specific antibody. The section was pre-treated using heat mediated antigen retrieval with Tris-EDTA buffer (pH 8.0-8.4) for 20 minut

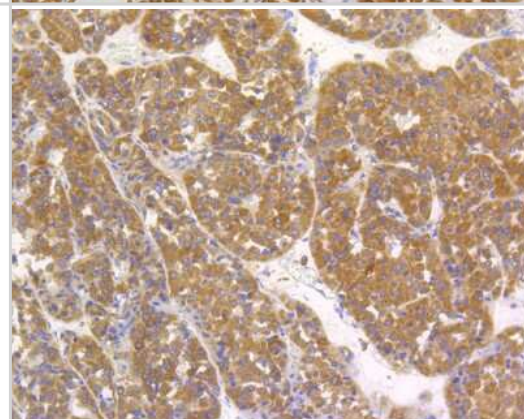

Immunohistochemistry-Paraffin: Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00) [NBP2-67295] - Immunohistochemical analysis of paraffin-embedded mouse kidney tissue using anti-Alkaline Phosphatase, Tissue Non-Specific antibody. The section was pre-treated using heat mediated antigen retrieval with Tris-EDTA buffer (pH 8.0-8.4) for 20 minutes. The t

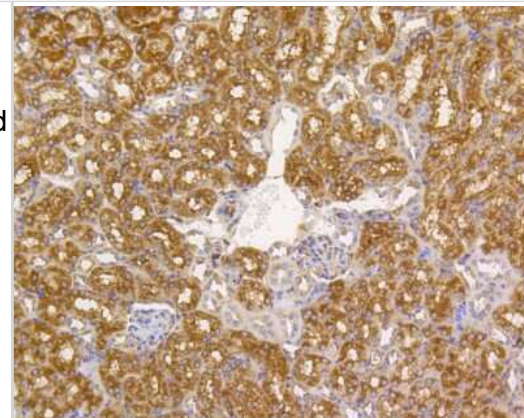

Western Blot: Alkaline Phosphatase, Tissue Non-Specific Antibody (SA40-00) [NBP2-67295] -Analysis of Alkaline Phosphatase on different lysates with Rabbit anti-Alkaline Phosphatase antibody at 1/5,000 dilution. Lane 1: Saos-2 cell lysate (15ug/Lane) Lane 2: HeLa cell lysate (15 ug/Lane) Lane 3: A549 cell lysate (15 ug/Lane) Lane 4: Mouse liver tissue lysate (20 ug/Lane) Lane 5: Rat liver tissue lysate (20 ug/Lane) Predicted band size: 57 kDa Observed band size: 75 kDa Exposure time: 24 seconds; 4-20% SDS-PAGE gel. Proteins were transferred to a PVDF membrane and blocked with 5% NFDm/TBST for 1 hour at room temperature. The primary antibody at 1/5,000 dilution was used in 5% NFDm/TBST at 4°C overnight. Goat Anti-Rabbit IgG - HRP Secondary Antibody at 1:50,000 dilution was used for 1 hour at room temperature.

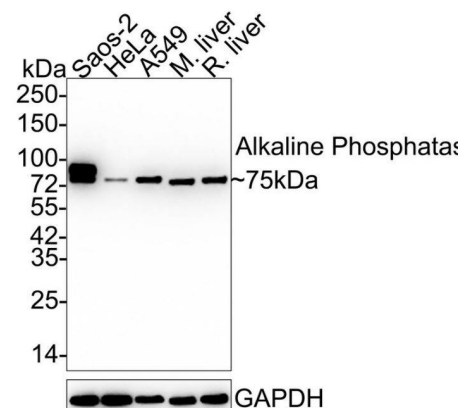

Publications

F Wei, CJ Neal, TS Sakthivel, Y Fu, M Omer, A Adhikary, S Ward, KM Ta, S Moxon, M Molinari, J Asiatico, M Kinzel, SN Yarmolenko, V San Cheong, N Orlovskaya, R Ghosh, S Seal, M Coathup A novel approach for the prevention of ionizing radiation-induced bone loss using a designer multifunctional cerium oxide nanozyme Bioactive materials, 2022-09-21;21(0):547-565. 2022-09-21 [PMID: 36185749]

Katsuhiro Kawaai, Yukiko Kuroda, Koichi Matsuo Calcified Cartilage-Guided Identification of Osteogenic Molecules and Geometries ACS Biomaterials Science & Engineering 2024-04-18 [PMID: 38634615]

Meng Z, Zhang S, Li W et al. Cationic proteins from eosinophils bind bone morphogenetic protein receptors promoting vascular calcification and atherogenesis European heart journal 2023-06-03 [PMID: 37279475] (WB, Mouse, Human)

Wang S, Gong Y, Wang Z et al. Regulon active landscape reveals cell development and functional state changes of human primary osteoblasts in vivo Human genomics 2023-02-15 [PMID: 36793138] (IHC-Fr)

Fan C, Li Z, Ji Q et al. Carboxymethyl chitin or chitosan for osteoinduction effect on the human periodontal ligament stem cells Dental materials journal 2022-03-04 [PMID: 35249903] (ICC/IF, Human)

Sun Y, Li R, Yu X Et Al. Highly active biological dermal acellular tissue scaffold composite with human bone powder for bone regeneration Materials & Design 2021-11-01 (ICC/IF)

Wei F, Neal C, Sakthivel T et al. Cerium oxide nanoparticles protect against irradiation-induced cellular damage while augmenting osteogenesis Materials Science and Engineering: C 2021-04-01 [PMID: 34082956]

Jiang L, Ding S, Ding W et al. Injectable sericin based nanocomposite hydrogel for multi-modal imaging-guided immunomodulatory bone regeneration Chemical Engineering Journal 2021-08-01 (WB, ICC/IF)



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Products Related to NBP2-67295

HAF008	Goat anti-Rabbit IgG Secondary Antibody [HRP]
NB7160	Goat anti-Rabbit IgG (H+L) Secondary Antibody [HRP]
NBP2-24891	Rabbit IgG Isotype Control
NBP1-91660PEP	Alkaline Phosphatase, Tissue Non-Specific Recombinant Protein Antigen

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-67295

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

