Product Datasheet

Colon Matched Tumor and Normal Tissue Slides (Tumor)Frozen NBP2-77670

Unit Size: 10 Slides

Store at -70C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP2-77670

Updated 10/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NBP2-77670



NBP2-77670


Colon Matched Tumor and Normal Tissue Slides (Tumor)- Frozen

Colon Matched Turnor and Normal Tissue Sildes (Turnor)- Frozen	
Product Information	
Unit Size	10 Slides
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -70C. Avoid freeze-thaw cycles.
Product Description	
Description	Frozen tissue slides are ideal for rapid cellular localization of RNA and protein expression, especially for antibodies that don't recognize paraffin embedded tissue. Tissues were snap frozen in liquid nitrogen immediately after excision, and embedded in OCT. Frozen tissue sections in each pair are prepared from the same donor. Tissue sections are 5-10 m in thickness, mounted on positively charged glass slides, and fixed by cold acetone. At least one of the tissue slides from each lot was stained with H & E for quality assurance. Product is supplied as 5 diseased slides and 5 adjacent normal slides.
Species	Human
Notes	Donor information available upon request
Lysate Type	Matched Tumor and Normal Tissue
Lysate Tissue	Colon
Lysate Tissue Condition	Tumor
Product Application Details	
Applications	Hematoxylin and Eosin Stain, Immunohistochemistry, Immunohistochemistry-Frozen, In-situ Hybridization
Recommended Dilutions	Immunohistochemistry, Immunohistochemistry-Frozen, In-situ Hybridization, Hematoxylin and Eosin Stain
Application Notes	Each individual tissue spot should be >90% core retention and total core lost from 1 slide should be no more than 5%
·	




Images

Hematoxylin & Eosin Stain: Colon Matched Tumor & Normal Tissue Slides (Tumor) [NBP2-77670] - Tissue: Human Colon Tumor, Pathology: Primary Tumor

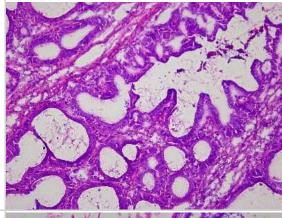

Hematoxylin & Eosin Stain: Colon Matched Tumor & Normal Tissue Slides (Tumor) [NBP2-77670] - Tissue: Human Colon, Pathology: Adjacent normal tissue

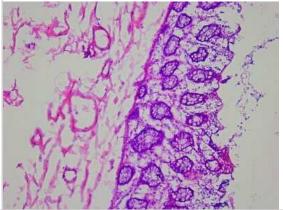



Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@bio-

techne.com
Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Products Related to NBP2-77670

NBP2-30256

Human Colon Tissue MicroArray (Rectal Cancer)

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Slides are guaranteed for 3 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-77670

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

