Product Datasheet

Kidney Tissue Slides (Papillary Carcinoma)- Paraffin NBP2-77905

Unit Size: 5 Slides

Store at 4C.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP2-77905

Updated 10/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NBP2-77905



NBP2-77905

Kidney Tissue Slides (Papillary Carcinoma)- Paraffin


Than by Thouse Chaos (Fapinary Caromona) Taramir		
Product Information		
Unit Size	5 Slides	
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.	
Storage	Store at 4C.	
Product Description		

Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at 4C.
Product Description	
Description	Each slide contains a single tissue section with 5 um thickness that is mounted on a positively charged glass slide. The slides included in this package are adjacent/serial sections of kidney tissue from one donor. Tissue was fixed in formalin immediately after excision and embedded in paraffin. This product can be used for both immunohistochemistry and in-situ hybridization. At least one of the tissue slides from each lot was stained with H&E to ensure the quality (not included in the package). Samples are IRB-approved from consented donors. Documentation on tissues' clinical histories may be available upon request. Donor information is also available upon request. Please contact nb-technical@bio-techne.com for any questions or requests. For FFPE DNA and RNA isolation, we may be able to provide tissue section slices in tubes. Please contact nb-custom@bio-techne.com for all custom requests related to this product.
Species	Human
Notes	Donor information available upon request
Lysate Type	Tissue
Lysate Tissue	Kidney
Lysate Tissue Condition	Panillany Carcinoma

Lysate Tissue Condition	Papillary Carcinoma
Product Application Details	
Applications	Immunohistochemistry, Immunohistochemistry-Paraffin, In-situ Hybridization, Dual RNAscope ISH-IHC
Recommended Dilutions	Immunohistochemistry, Immunohistochemistry-Paraffin, In-situ Hybridization, Dual RNAscope ISH-IHC
Application Notes	Please bake slides at 60C for 30 minutes before use.

Images

Hematoxylin & Eosin Stain: Kidney Tissue Slides (Papillary Carcinoma) [NBP2-77905] - Tissue: Human Kidney, Pathology: Papillary Renal Cell Carcinoma

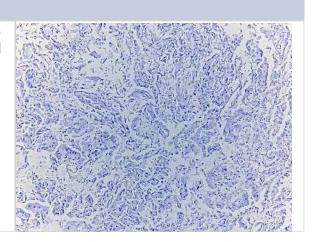





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Products Related to NBP2-77905

NBP2-30220 Human Kidney Tissue MicroArray (Cancer)

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Slides are guaranteed for 3 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-77905

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

