

Product Datasheet

Exosome Standards (A549 cell line) NBP3-11645-180ug

Unit Size: 6 x 30ug Vials

Store at 4C. After reconstitution store at -70C.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP3-11645

Updated 4/13/2025 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP3-11645



NBP3-11645-180ug

Exosome Standards (A549 cell line)

Product Information	
Unit Size	6 x 30ug Vials
Concentration	Please see the protocols for proper use of this product. If no protocol is available, contact technical services for assistance.
Storage	Store at 4C. After reconstitution store at -70C.
Reconstitution Instructions	Add deionized water, 100 ul for Standard 100 ug and 30 ul for Standard 30 ug, to get a final concentration of 1 mg/mL. Resuspend exosomes pipetting the solution up and down 10-15 times, avoiding bubbles. Vortex the reconstituted standard for 60 seconds.
Buffer	Lyophilized from cell culture media


Product Description	
Description	Highly pure, lyophilized exosome standards with superior stability, optimal for multiple applications including: • Assay calibration • Spike-in control for exosome quantification • Protein marker analysis for different techniques such as Western Blot and Flow Cytometry • Extraction and analysis of exosomal RNA and DNA • Standardized positive controls for evaluating immunocapture performance Quantity per vial of 30 ug size (number of particles in 30 ug: $> 1 \times 10^8$). Quantity per vial of 100 ug size (number of particles in 100 ug: $> 1 \times 10^{10}$).
Preparation Method	Isolation involves Tangential flow filtration combined with Size Exclusion Chromatography. Exosomes (small EVs) are quantified and validated for protein content and particle number by Nanoparticle Tracking Analysis as well as for common tetraspanins marker validation. Lyophilization does not alter stability of exosome proteins and nucleic acids.

Product Application Details	
Applications	Simple Western, ELISA, Electron Microscopy, Flow Cytometry, Nucleic Acid Extraction
Recommended Dilutions	Simple Western, Flow Cytometry, ELISA, Electron Microscopy, Nucleic Acid Extraction




Images


Flow Cytometry: Exosome Standards (A549 cell line) [NBP3-11645] - Phenotyping assays by FACS. Reconstituted Exosomes can be used for profiling biomarkers by FACS analysis. Recommended quantity: 5 ug of reconstituted Exosomes Standards for each test.

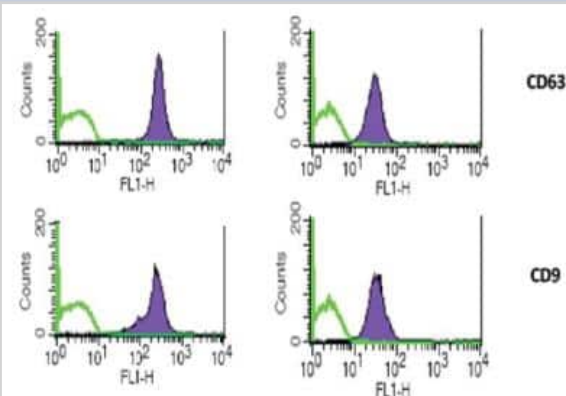

Simple Western: Exosome Standards (A549 cell line) [NBP3-11645] - Simple Western lane view shows lysates of Exosome Standards (A549 cell line), exosome standards (HT29 cell line), and COLO 205 human colorectal adenocarcinoma cell line whole cell lysate (WCL), loaded at 0.2 mg/mL. A specific band was detected for Annexin V at approximately 38 kDa (as indicated) using 50 ug/mL of Mouse Anti-Human Annexin V Monoclonal Antibody. This experiment was conducted under reducing conditions and using the 12-230 kDa separation system.

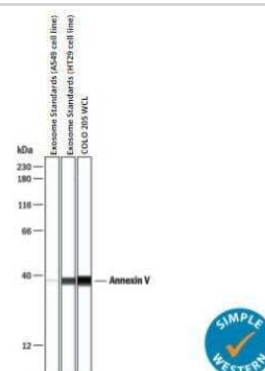

Simple Western: Exosome Standards (A549 cell line) [NBP3-11645] - Simple Western lane view shows lysates of Exosome Standards (A549 cell line), exosome standards (HT29 cell line), and COLO 205 human colorectal adenocarcinoma cell line whole cell lysate (WCL), loaded at 0.2 mg/mL. A specific band was detected for Flotillin-1 at approximately 51 kDa (as indicated) using 1:250 ug/mL of Flotillin-1 Antibody (NBP1-79022). This experiment was conducted under reducing conditions and using the 12-230 kDa separation system.

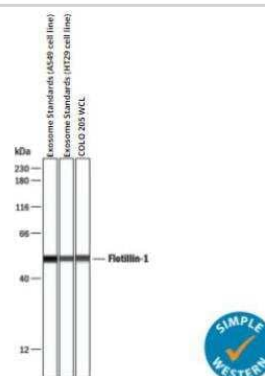



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Support products are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-11645

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

