Product Datasheet

Exosome Standards (LnCAP cell line) NBP3-11687-200ug

Unit Size: 2 x 100ug Vials

Store at 4C. After reconstitution store at -70C.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP3-11687

Updated 4/2/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NBP3-11687



NBP3-11687-200ug


Exosome Standards (LnCAP cell line)

Exosome Standards (Enorth Son line)	
Product Information	
Unit Size	2 x 100ug Vials
Concentration	Please see the protocols for proper use of this product. If no protocol is available, contact technical services for assistance.
Storage	Store at 4C. After reconstitution store at -70C.
Reconstitution Instructions	Add deionized water, 100 ul for Standard 100 ug and 30 ul for Standard 30 ug, to get a final concentration of 1 mg/mL. Resuspend exosomes pipetting the solution up and down 10-15 times, avoiding bubbles. Vortex the reconstituted standard for 60 seconds.
Buffer	Lyophilized from cell culture media
Product Description	
Description	Highly pure, lyophilized exosome standards with superior stability, optimal for multiple applications including: Assay calibration, Spike-in control for exosome quantification, Protein marker analysis for different techniques such as Western Blot and Flow Cytometry, Extraction and analysis of exosomal RNA and DNA, Standardized positive controls for evaluating immunocapture performance. Quantity per vial of 30 ug size (number of particles in 30 ug: > 1x10^8). Quantity per vial of 100 ug size (number of particles in 100 ug: > 1x10^10). Store lyophilized exosomes up to 3 years at 4C (Expire date indicated on the label). For reconstituted exosomes store at -20C for up to one month or at -80C for up to 6 months. Avoid repeated freeze-and-thraw cycles.
Preparation Method	Isolation involves Tangential flow filtration combined with Size Exclusion Chromatogrpahy. Exosomes (small EVs) are quantified and validated for protein content and particle number by Nanoparticle Tracking Analysis as well as for common tetraspanins maker validation. Lyophilization does not alter stability of exosome proteins and nucleic acids.

Product Application Details	
Applications	ELISA, Electron Microscopy, Flow Cytometry, Nucleic Acid Extraction
Recommended Dilutions	Flow Cytometry, ELISA, Electron Microscopy, Nucleic Acid Extraction

Images

Flow Cytometry: Exosome Standards (LnCAP cell line) [NBP3-11687] - Phenotyping assays by FACS. Reconstituted Exosomes can be used for profiling biomarkers by FACS analysis. Recommended quantity: 5 ug of reconstituted Exosomes Standards for each test.

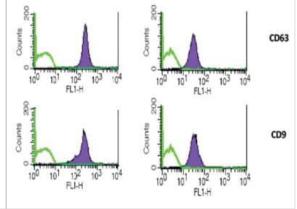



Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Support products are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-11687

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

