Product Datasheet

TIMP-2 Antibody (TIMP2/8193R) NBP3-20229-100ug

Unit Size: 100 ug

Store at 4C. Do not freeze.

www.novusbio.com

technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP3-20229

Updated 7/16/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NBP3-20229



NBP3-20229-100ug


TIMP-2 Antibody (TIMP2/8193R)

Product Information	
Unit Size	100 ug
Concentration	0.2 mg/ml
Storage	Store at 4C. Do not freeze.
Clonality	Monoclonal
Clone	TIMP2/8193R
Preservative	0.05% Sodium Azide
Isotype	IgG Kappa
Purity	Protein A or G purified
Buffer	10mM PBS with 0.05% BSA
Product Description	
Description	Positive Controls: Normal Placenta, Prostate or Ovarian Carcinoma.
	Antibody with azide - store at 2 to 8C. Antibody without azide - store at -20 to -80C. Non-hazardous. No MSDS required.
Host	Rabbit
Gene ID	7077
Gene Symbol	TIMP2
Species	Human
Specificity/Sensitivity	It recognizes a protein of 21kDa, identified as tissue inhibitor of metalloproteinases-2 (TIMP-2) It is closely related to TIMP-1 and shows the highest binding affinity to both the latent (pro) and active forms of 72kDa Type IV collagenase (also known as MMP-2 or gelatinase A) It also has affinity for the active form of 92kDa Type IV collagenase (also known as MMP-9 or gelatinase B) TIMP-2 inhibits the proteolytic invasiveness of tumor cells and normal placental trophoblast cells
Immunogen	Recombinant fragment of human TIMP-2 protein (exact sequence is proprietary)
Product Application Details	
Applications	Immunohistochemistry-Paraffin
Recommended Dilutions	Immunohistochemistry-Paraffin 1-2 ug/ml
Application Notes	Immunohistochemistry (Formalin-fixed): 1-2ug/ml for 30 minutes. at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95C followed by cooling at RT for 20 minutes.
	oplimation for a specific application should be determined.



Images

IHC analysis of formalin-fixed, paraffin-embedded human brain. Negative tissue control using TIMP-2 antibody (TIMP2/8193R) at 2ug/ml in PBS for 30min RT. HIER: Tris/EDTA, pH9.0, 45min. Secondary: HRP-polymer, 30min. DAB, 5min.

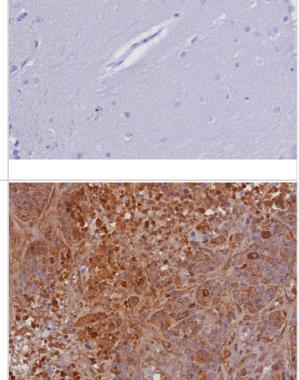

Formalin-fixed, paraffin-embedded human bladder carcinoma stained with TIMP-2 antibody (TIMP2/8193R). HIER: Tris/EDTA, pH9.0, 45min. Secondary: HRP-polymer, 30min. DAB, 5min.



www.novusbio.com



Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA Phone: 303.730.1950 Toll Free: 1.888.506.6887 Fax: 303.730.1966 nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402 canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449 Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com Orders: nb-customerservice@bio-techne.com General: novus@novusbio.com

Products Related to NBP3-20229-100ug

210-TA-005	TNF-alpha [Unconjugated]
NBP2-61356-10ug	Recombinant Human TIMP-2 Protein
NB7160	Goat anti-Rabbit IgG (H+L) Secondary Antibody [HRP]
HAF008	Goat anti-Rabbit IgG Secondary Antibody [HRP]

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-20229

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

www.novusbio.com

