# **Product Datasheet**

# Fibronectin Antibody (CL3730) - Azide and BSA Free NBP3-43942

Unit Size: 100 ul

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP3-43942

Updated 3/13/2025 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NBP3-43942



# NBP3-43942

Fibronectin Antibody (CL3730) - Azide and BSA Free

Tibionectin Antibody (CE3730) - Azide and BSA Tiee	
Product Information	
100 ul	
LYOPH mg/ml	
Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.	
Monoclonal	
CL3730	
No Preservative	
Centrifuge the vial of lyophilized antibody at 12,000 x g for 20 seconds. Reconstitute by adding sterile, distilled water to achieve a final antibody concentration of 1mg/ml.	
IgG1	
Protein A purified	
Lyophilized from a 0.2 um filtered solution in PBS with Trehalose	
Mouse	
2335	
FN1	
Human	
This antibody was generated using a recombinant protein sequence of P02751, with the exact immunogen sequence remaining proprietary.	
Product Application Details	
Western Blot, Immunohistochemistry-Paraffin	
Western Blot 1 ug/ml, Immunohistochemistry-Paraffin 1:1000 - 1:2500	


IHC-Paraffin HIER pH6 retrieval is recommended.

# **Images**

**Application Notes** 

Lane 1: Marker [kDa] 250, 130, 100, 70, 55, 35, 25, 15, 10

Lane 2: Human cell line U-87 MG

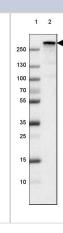

Page 2 of 4 v.20.1 Updated 3/13/2025 Staining of human breast cancer shows strong immunoreactivity in tumor stroma. Staining of human liver cancer shows strong positivity in tumor. Staining of human endometrium shows strong immunoreactivity in endometrial stroma and basement membrane of glandular cells. Staining of human prostate shows moderate positivity in smooth muscle.



Page 3 of 4 v.20.1 Updated 3/13/2025 Staining of human placenta shows strong immunoreactivity in chorionic Staining of human pancreas shows absence of immunoreactivity in glandular cells (negative control). Staining of human testis shows immunoreactivity in the connective tissue surrounding seminiferous tubules. Staining of human small intestine shows strong immunoreactivity in the connective tissue.





### **Novus Biologicals USA**

10730 E. Briarwood Avenue Centennial, CO 80112

USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

#### **Bio-Techne Canada**

21 Canmotor Ave Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

#### **Bio-Techne Ltd**

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

#### **General Contact Information**

www.novusbio.com

Technical Support: nb-technical@bio-

techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

# **Products Related to NBP3-43942**

HAF007 Goat anti-Mouse IgG Secondary Antibody [HRP]

NB720-B Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin]

NBP1-97005-0.5mg Mouse IgG1 Isotype Control (MG1)
NBP1-91258PEP Fibronectin Antibody Blocking Peptide

#### Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-43942

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

