Product Datasheet

NECAB1 Antibody (CL0576) - Azide and BSA Free NBP3-43976

Unit Size: 100 ul

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP3-43976

Updated 3/11/2025 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NBP3-43976




NBP3-43976

NECAB1 Antibody (CL0576) - Azide and BSA Free


THEORET Filliability (CESOTO) Filliability 100	
Product Information	
Unit Size	100 ul
Concentration	LYOPH mg/ml
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	CL0576
Reconstitution Instructions	Centrifuge the vial of lyophilized antibody at 12,000 x g for 20 seconds. Reconstitute by adding sterile, distilled water to achieve a final antibody concentration of 1mg/ml.
Isotype	lgG2b
Purity	Protein A purified
Buffer	Lyophilized from a 0.2 um filtered solution in PBS with Trehalose
Product Description	
Host	Mouse

-	
Host	Mouse
Gene ID	64168
Gene Symbol	NECAB1
Species	Human
Immunogen	This antibody was generated using a recombinant protein sequence of Q8N987, with the exact immunogen sequence remaining proprietary.
Product Application Dataila	


Product Application Details	
Applications	Western Blot, Immunohistochemistry-Paraffin
Recommended Dilutions	Western Blot 1 ug/ml, Immunohistochemistry-Paraffin 1:500 - 1:1000
Application Notes	For IHC-Paraffin, HIER pH 6 retrieval is recommended.

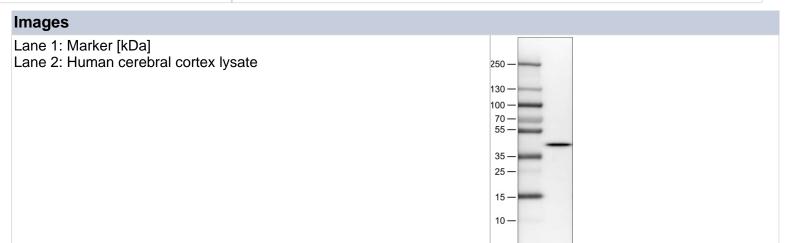


Staining of human cerebellum shows moderate cytoplasmic positivity in neurons.

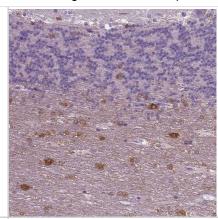


Analysis in human cerebral cortex and pancreas tissues using NBP3-43976 antibody. Corresponding NECAB1 RNA-seq data are presented for the same tissues.

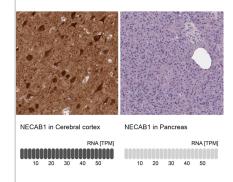


Staining of human cerebral cortex shows strong cytoplasmic positivity in neurons.

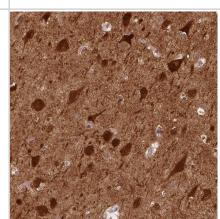

Staining of human kidney shows strong cytoplasmic positivity in cells in glomeruli.

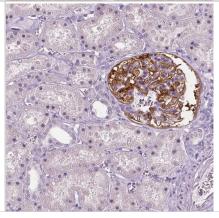

Staining of human pancreas shows no positivity in exocrine glandular cells as expected.

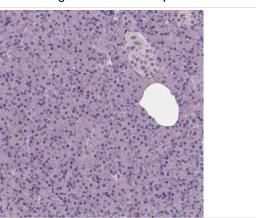



Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-43976

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

