# **Product Datasheet**

# Histamine H3R Antibody - BSA Free NLS3262

Unit Size: 0.05 ml

Keep as concentrated solution. Aliquot and store at -20C or below. Avoid multiple freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NLS3262

Updated 10/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NLS3262



# **NLS3262**


Histamine H3R Antibody - BSA Free

nistamine har Antibody - Bar Free	
Product Information	
0.05 ml	
1.0 mg/ml	
Keep as concentrated solution. Aliquot and store at -20C or below. Avoid multiple freeze-thaw cycles.	
Polyclonal	
0.1% Sodium Azide	
IgG	
Immunogen affinity purified	
PBS	
Product Description	
Product can be stored undiluted at 4C for up to 1 month.	
Rabbit	
11255	
HRH3	
Human, Mouse, Rat, Bovine, Equine, Guinea Pig, Hamster, Primate, Monkey	
Predicted cross-reactivity based on sequence identity: Gorilla (100%), Gibbon (100%), Opossum (81%).	
Human Histamine H3R. BLAST analysis of the peptide immunogen showed no homology with other human proteins.	
Synthetic 16 amino acid peptide from 2nd cytoplasmic domain of human Histamine H3R.	
Product Application Details	
Immunohistochemistry, Immunohistochemistry-Paraffin, Western Blot (Negative)	
Immunohistochemistry, Immunohistochemistry-Paraffin 7 - 15 ug/ml, Western Blot (Negative)	




# **Images**

Immunohistochemistry-Paraffin: Histamine H3 R Antibody [NLS3262] - Analysis of anti-Histamine 3 Receptor / HRH3 antibody with human brain, neurons and glia at dilution 7-15 ug/ml.

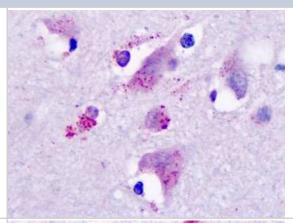

Immunohistochemistry-Paraffin: Histamine H3 R Antibody [NLS3262] - Analysis of anti-Histamine 3 Receptor / HRH3 antibody with human brain, neurons.

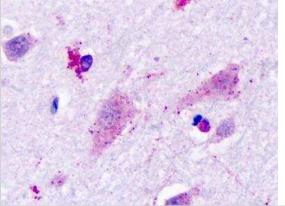

#### **Procedures**

#### Immunohistochemistry Protocol for Histamine receptor H3 Antibody (NLS3262)

Immunohistochemistry Protocol for Histamine receptor H3 Antibody (NLS3262): Immunohistochemistry

- 1. Prepare tissue with formalin fixation and by embedding it in paraffin wax.
- 2. Make 4 um sections and place on pre-cleaned and charged microscope slides.
- 3. Heat in a tissue-drying oven for 45 minutes @ 60 degrees Celcius.
- 4. Deparaffinize the tissues by wash drying the slides in 3 changes of xylene for 5 minutes each @ RT.
- 5. Rehydrate the tissues by washing the slides in 3 changes of 100% alcohol for 3 minutes each @ RT.
- 6. Wash the slides in 2 changes of 95% alcohol for 3 minutes each @ RT.
- 7. Wash the slides in 1 change of 80% alcohol for 3 minutes @ RT.
- 8. Rinse the slides in gentle running distilled water for 5 minutes @ RT.
- Perform antigen retrieval by steaming the slides in 0.01M sodium citrate buffer (pH 6.0) @ 99-100 degrees Celcius for 20 minutes.
- 10. Remove the slides from the heat and let stand in buffer @ RT for 20 minutes.
- 11. Rinse the slides in 1X TBS-T for 1 minute @ RT.
- \*\*Do not allow the tissues to dry at any time during the staining procedure\*\*
- 12. Begin the immunostaining by applying a universal protein block for 20 minutes @ RT.
- 13. Drain protein block from the slides and apply the diluted primary antibody for 45 minutes @ RT.
- 14. Rinse the slide in 1X TBS-T for 1 minute @ RT.
- 15. Apply a biotinylated anti-rabbit IgG (H+L) secondary for 30 minutes @ RT.
- 16. Rinse the slide in 1X TBS-T for 1 minute @ RT.
- 17. Apply an alkaline phosphatase steptavidin for 30 minutes @ RT.
- 18. Rinse the slide in 1X TBS-T for 1 minute @ RT.
- 19. Apply an alkaline phosphatase chromagen substrate for 30 minutes @ RT.
- 20. Rinse the slide in distilled water for 1 minute @ RT.
- \*\*This method should only be used if the chromagen substrate is alcohol insoluble (ie: Vector Red, DAB)\*\*
- 21. Dehydrate the tissue by washing the slides in 2 changes of 80% alcohol for 1 minute each @ RT.
- 22. Wash the slides in 2 changes of 95% alcohol for 1 minute each @ RT.
- 23. Wash the slides in 3 changes of 100% alcohol for 1 minute each @ RT.
- 24. Wash the slides in 3 changes of xylene for 1 minute each @ RT.
- 25. Apply cover slip.





# **Novus Biologicals USA**

10730 E. Briarwood Avenue Centennial, CO 80112

USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449 Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

## **General Contact Information**

www.novusbio.com

Technical Support: nb-technical@bio-

techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

# **Products Related to NLS3262**

HAF008 Goat anti-Rabbit IgG Secondary Antibody [HRP]

NB7160 Goat anti-Rabbit IgG (H+L) Secondary Antibody [HRP]

NBP2-24891 Rabbit IgG Isotype Control 210-TA-005 TNF-alpha [Unconjugated]

#### Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NLS3262

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

